


Antibacterials

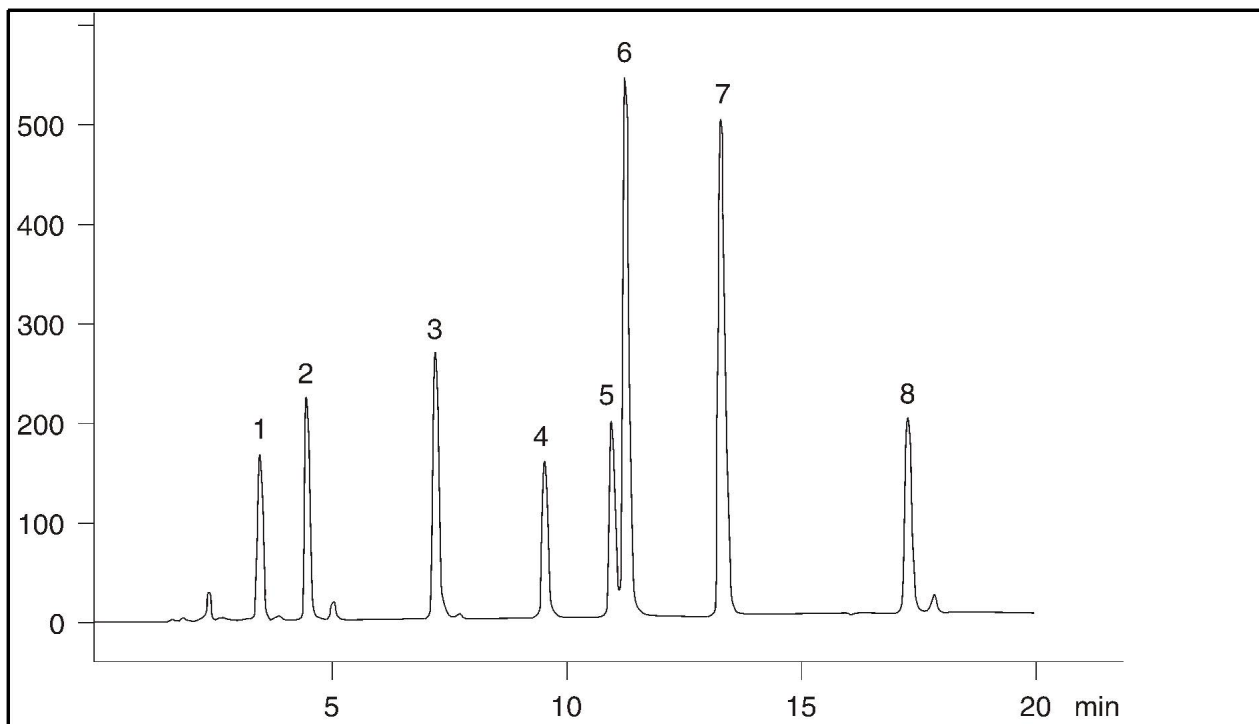

Column: Luna 5 μ m C6H5 / Phenyl Hexyl

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Gradient

Eluent A: 20 mM KH₂PO₄

Eluent B: Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	80	20
2	5	75	25
3	20	55	45

Flow Rate: 1 mL/min

Temperature: 40 °C

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm

SPE:

SAMPLE ANALYTE

- 1 Carbadox
- 2 Thiamphenicol
- 3 Furazolidone
- 4 Oxolinic acid
- 5 Sulfadimethoxine
- 6 Sulfaquinoxaline
- 7 Nalidixic acid

8 Piromiodic acid

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1201/ALG/50