


APPLICATION NOTE - HPLC
Tricyclic Antidepressants at pH7.0



15029

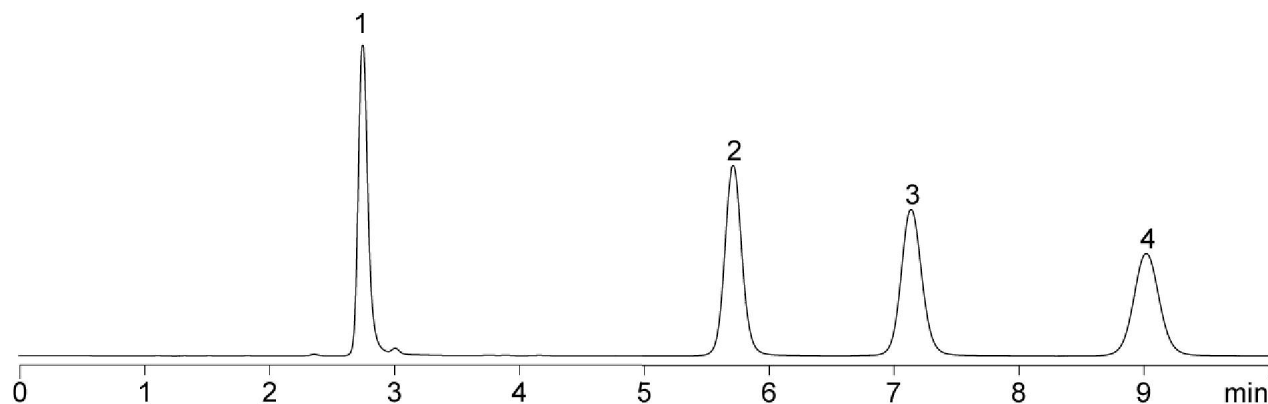

Column: Gemini 5 μ m C18 / ODS C18

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4435-E0

Elution Type: Isocratic

Eluent A: 20mM Phosphate Buffer, pH7.0

Pct. A: 30

Eluent B: ACN

Pct. B: 35

Eluent C: MeOH

Pct. C: 35

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTES:

- 1 Imipramine
- 2 Nortriptyline
- 3 Amitriptyline
- 4 Clomipramine

SPE:

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15029

Page 1 of 1

