


Benzodiazepines Standards on Gemini C18



15456

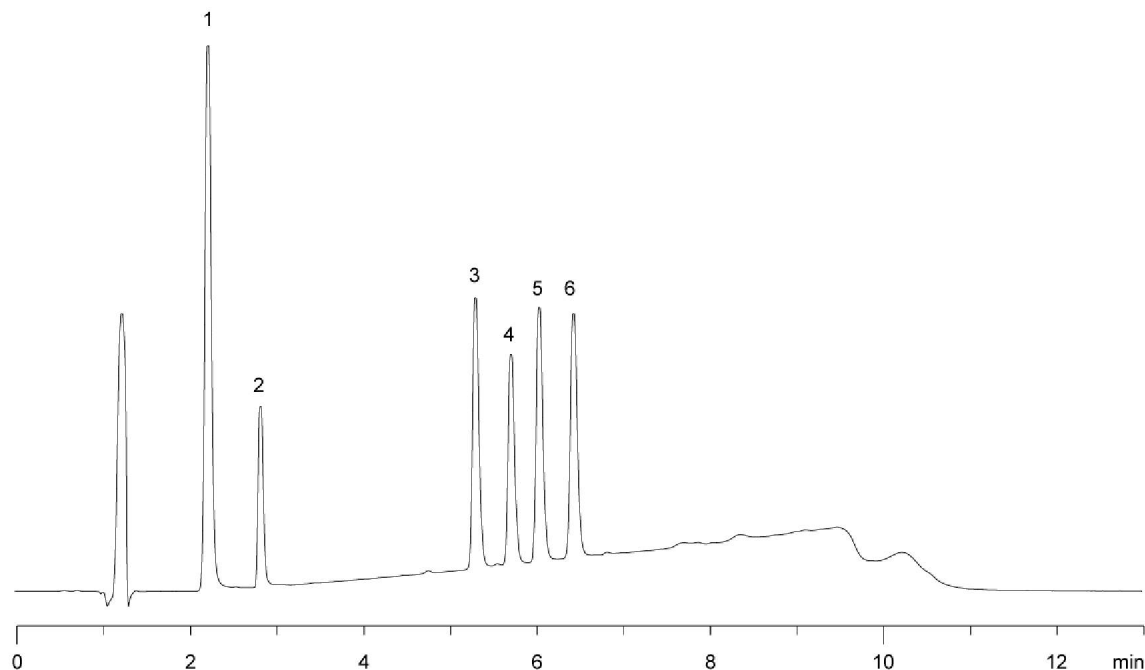

Column: Gemini 5 μ m C18 / ODS C18
Dimensions: 150 x 3 mm ID
Order No: 00F-4435-Y0
Elution Type: Gradient
Eluent A: 0.1% Formic acid in Water
Eluent B: 0.1% Formic acid in Acetonitrile

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	80	20
2	8	20	80
3	9	20	80

Flow Rate: 0.8 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm
SPE:

SAMPLE ANALYTES:

- 1 Chlordiazepoxide
- 2 Flurazepam
- 3 Oxazepam
- 4 Clonazepam
- 5 Temazepam
- 6 Diazepam

For more information call your Phenomenex Representative at

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 15456

Page 1 of 1