


Penicillin Antibiotics

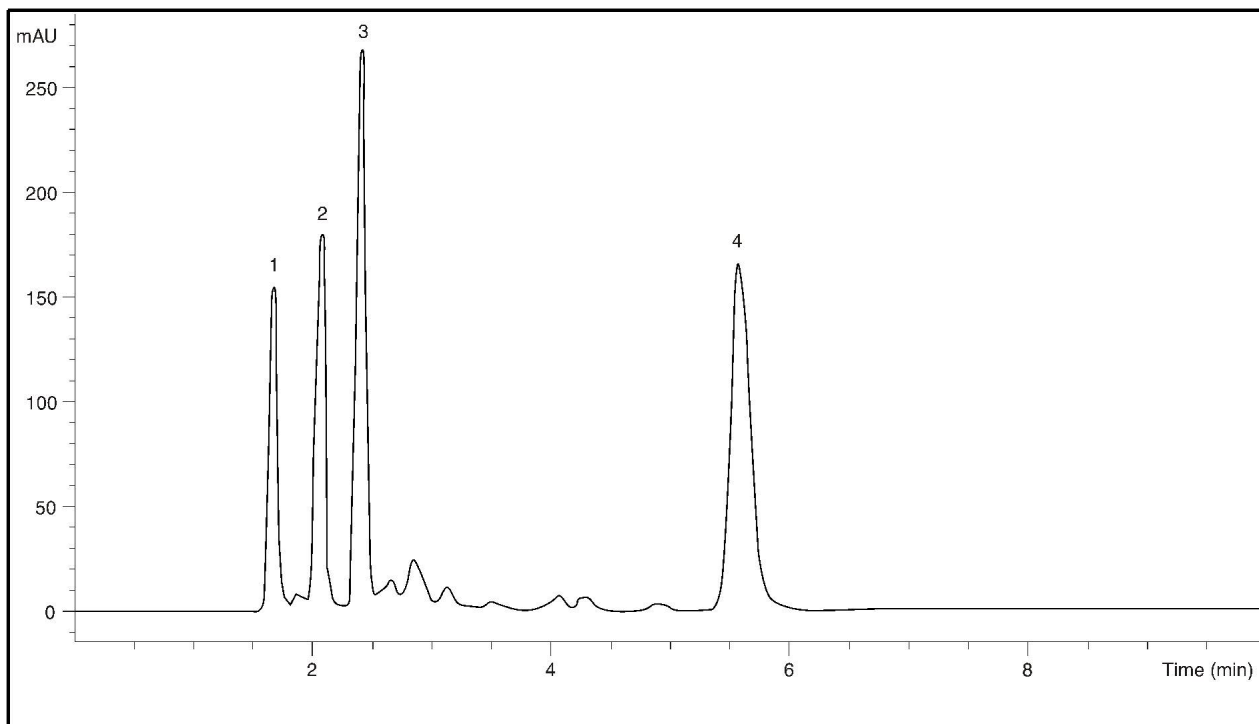

Column: Luna 5 μ m C8 Original

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 55

Eluent B: 25 mM KH₂PO₄, pH 4.6

Pct. B: 45

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 Amoxicillin
- 2 Ampicillin
- 3 Penicillin G
- 4 Penicillin V

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1215/ALG/68