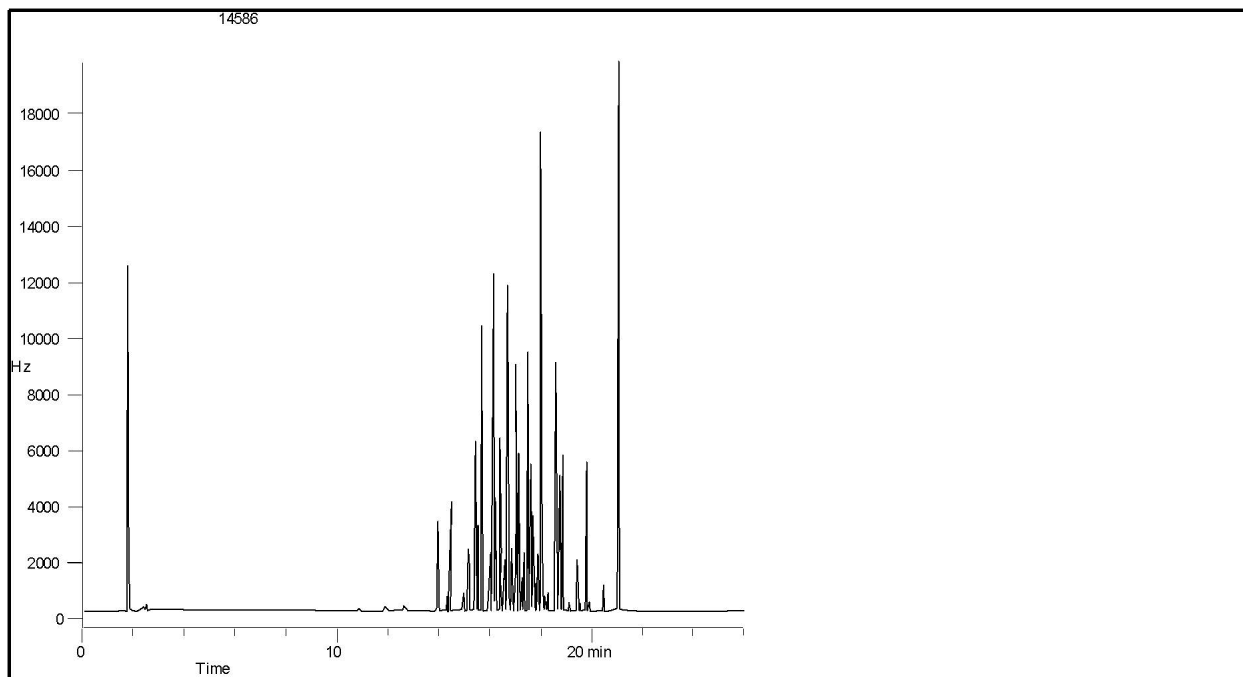


Aroclor 1260 on ZB-5 30m x 0.32mm x 0.25µm

Column: Zebron ZB5 30m x .32mm x .25um
Phase: 5%-diphenyl-95%-dimethyl polysiloxane copolymer
Dimensions: 30 meters x 0.32 mm x 0.25 µm
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	1.7	120
2	1	1.7	120
3	21	1.7	300
4	26	1.7	300

Carrier Gas: Helium 35.23 cm/sec
Flow Rate: 1.7 mL/min
Injection: Split 15:1 1 µL @ 250 °C
Detection: Electron Capture (ECD) 375 °C

SAMPLE ANALYTE

1 Aroclor 1260

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 14586