


EZ:faast Amino Acid Hydrolysate Standard (GC/MS)

Page 1 of 1


Column: EZ:faast GC
Phase: Zebron ZB-AAA (proprietary)
Dimensions: 10 meters x 0.25 mm x 0 µm
Order No: KG0-7168
Oven Profile: 30°C/min from 110° to 320°C, hold at 320°C for 1 minute
Carrier Gas: Constant Flow Helium, 1.1 mL/min
Injection: Split 15:1 2.5 µL @ 250 °C
Detection: Mass Spectrometer (MS) (300 °C)

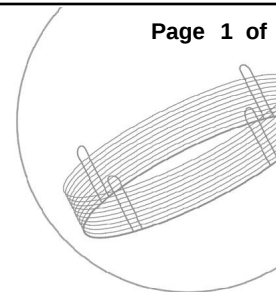

Products used in this Application:

EZ:faast™

14612

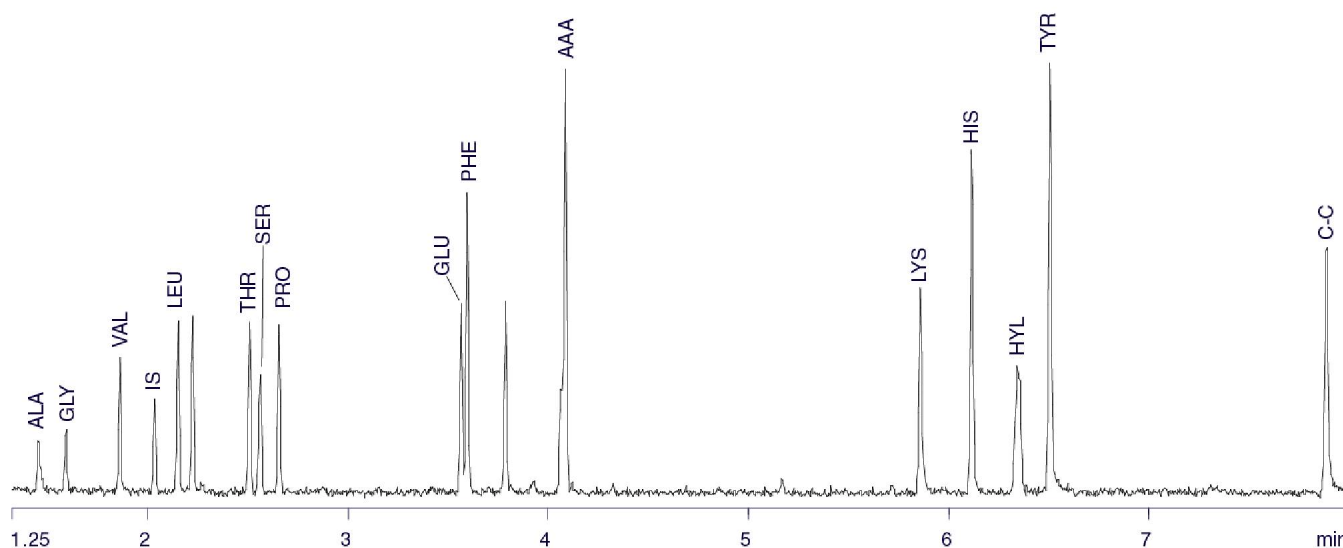

ANALYTES:

- | | |
|-------|--------|
| 1 ALA | 10 ASP |
| 2 GLY | 11 MET |
| 3 VAL | 12 HYP |
| 4 IS | 14 PHE |
| 5 LEU | 15 LYS |
| 6 ILE | 16 HIS |
| 7 THR | 17 HYL |
| 8 SER | 18 TYR |
| 9 PRO | 19 C-C |

©Copyright 2000-2005 Phenomenex Inc.

For more information call your Phenomenex Representative at (310) 212-0555



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


 ...breaking with tradition™

USA
 tel.: (310) 212-0555
 fax: (310) 328-7768
 email: info@phenomenex.com

Puerto Rico
 (800) 541-HPLC
 (310) 328-7768
 info@phenomenex.com

Canada
 (800) 543-3681
 (310) 328-7768
 info@phenomenex.com

France
 01 30 09 21 10
 01 30 09 21 11
 franceinfo@phenomenex.com

United Kingdom
 01625-501367
 01625-501796
 ukinfo@phenomenex.com

Germany
 06021-58830-0
 06021-58830-11
 anfrage@phenomenex.com

New Zealand
 09-4780951
 09-4780952
 info@phenomenex.co.nz

Australia
 02-9428-6444
 02-9428-6445
 info@phenomenex.com.au