


Antibiotics

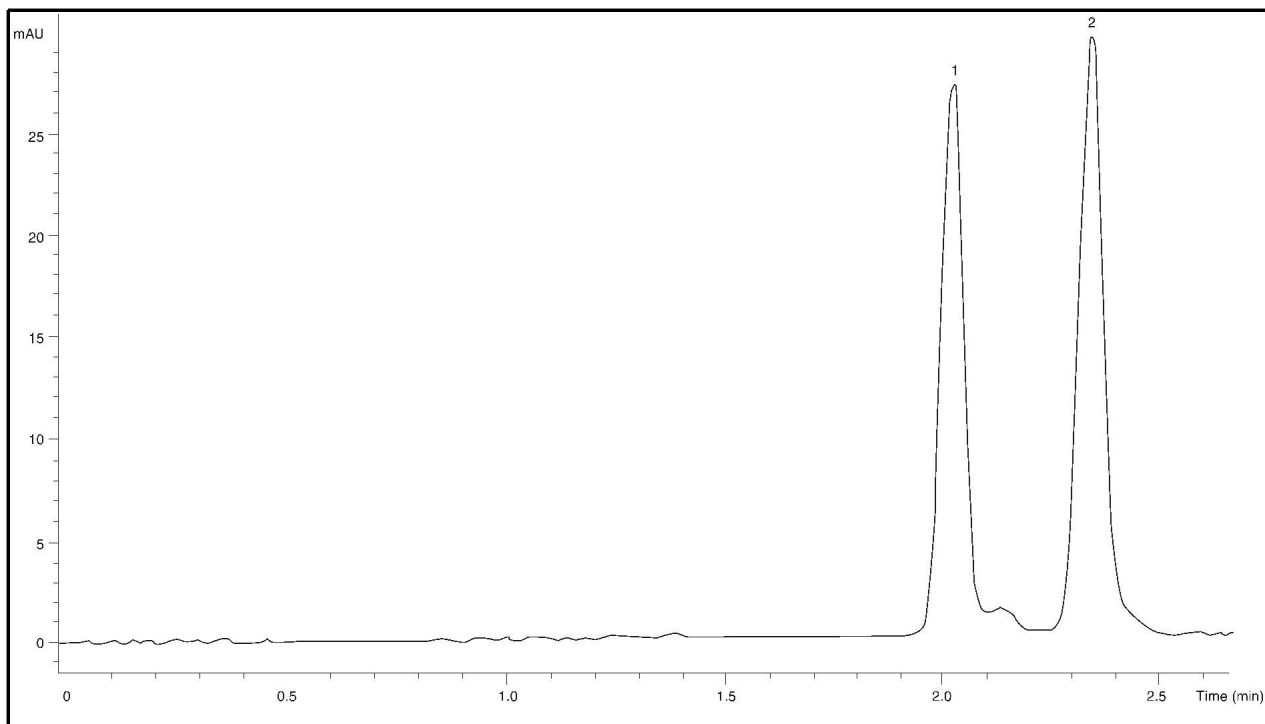

Column: Luna 5 μ m C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 30

Eluent B: Water

Pct. B: 70

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 365 nm

SAMPLE ANALYTE

1 Oxytetracycline

2 Tetracycline

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1245/ALG/106