

20mer Hetero-Oligonucleotide on Jupiter Proteo - ZOOM View

Page 1 of 1

Column: Jupiter 4u Proteo 90A

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4396-E0

Elution Type: Gradient

Eluent A: 0.1M TEAA, pH 7.0

Eluent B: Acetonitrile


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	20	85	15

Flow Rate: 1 mL/min

Col. Temp.: 60 °C

Detection: UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)

Products used in this Application:

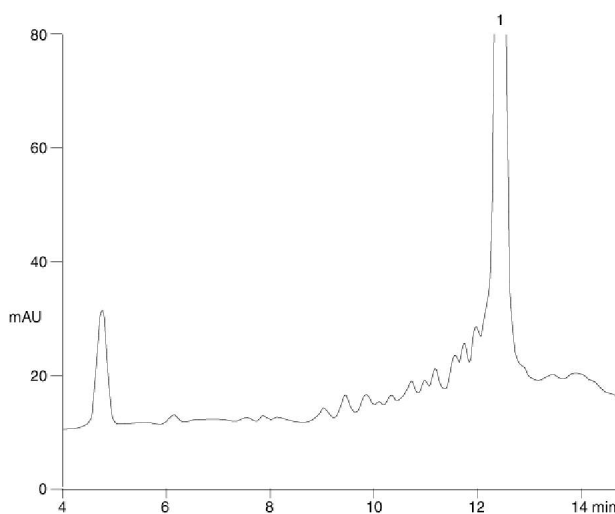

ANALYTES:

- 1 20mer heteroligonucleotide



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.