


Antibiotics

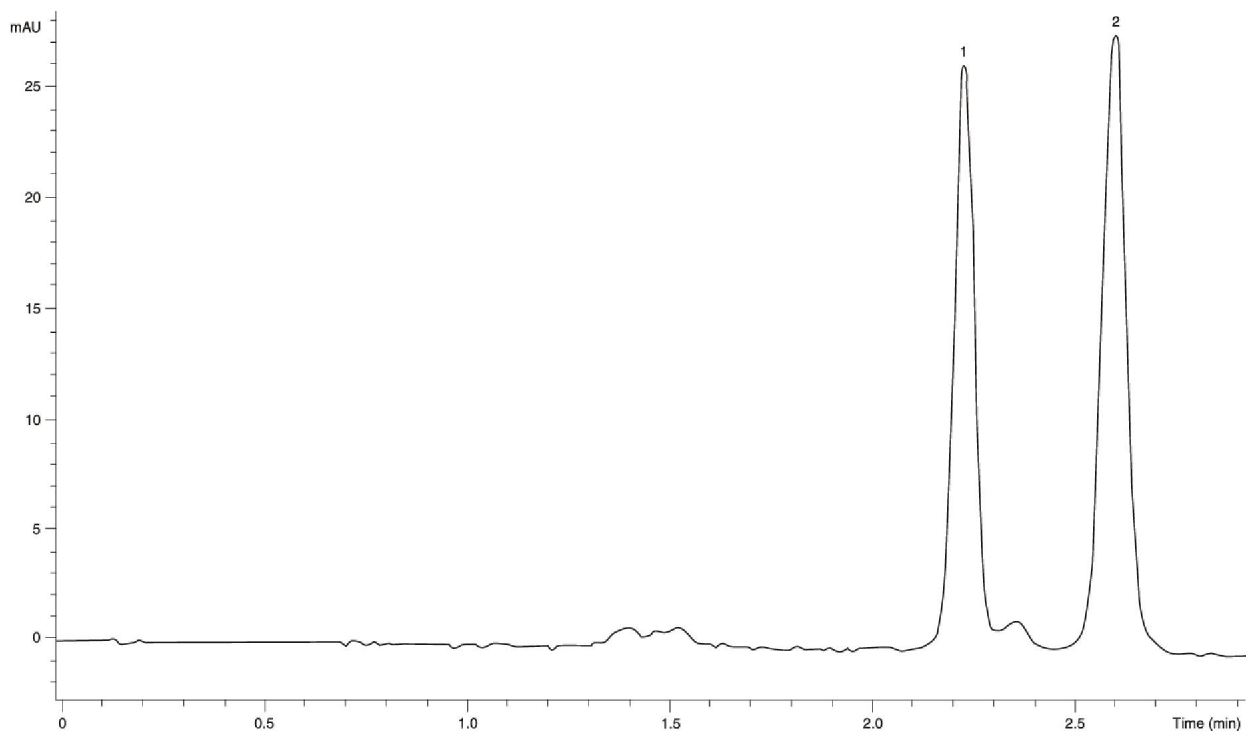

Column: Luna 5 μ m C8 (2)

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4249-E0

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 30

Eluent B: Water

Pct. B: 70

Flow Rate: 1 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 365 nm

SAMPLE ANALYTES:

1 Oxytetracycline

2 Tetracycline

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2005 Phenomenex Inc.

Ref. No.: 1246/ALG/107

Page 1 of 1