

## Herbicides

**Column:** Luna 3u C18(2) 100A

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4251-E0

**Elution Type:** Gradient

**Eluent A:** Acetonitrile


**Eluent B:** Water

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 70    | 30    |
|                   | 2        | 4          | 100   | 0     |

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

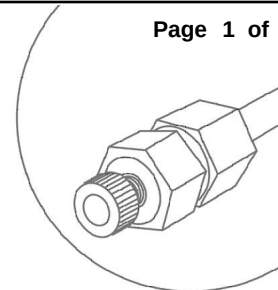

Products used in this Application:

**LUNA™**

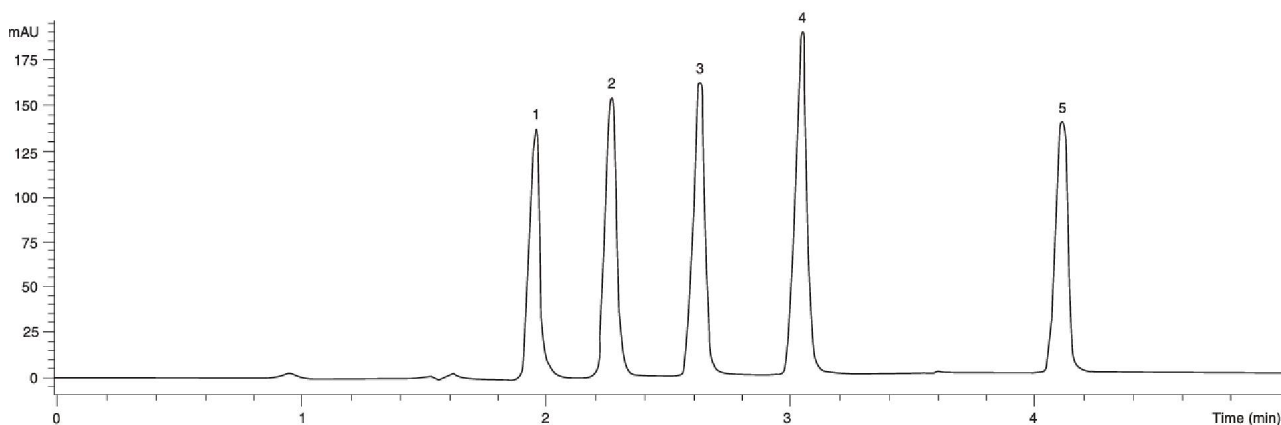

## ANALYTES:

- 1 Tebuthiron
- 2 Sinazine
- 3 Altrazine
- 4 Propazine
- 5 Dacthol

