


Ochratoxin analysis using Gemini

Page 1 of 1

Column: Gemini 5u C18 110A
Dimensions: 150 x 3 mm ID
Order No: 00F-4435-Y0
Elution Type: Isocratic
Eluent A: 10 mM Ammonium Acetate, pH 4.5
Eluent B: Acetonitrile
Flow Rate: 0.5 mL/min
Col. Temp.: ambient
Detection: Mass Spectrometer (MS) (ambient)

Pct. A: 50


Pct. B: 50

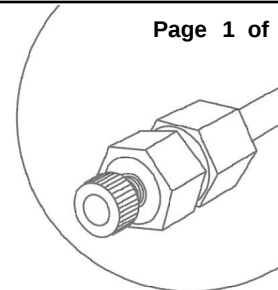

Products used in this Application:

Gemini[™]

15529

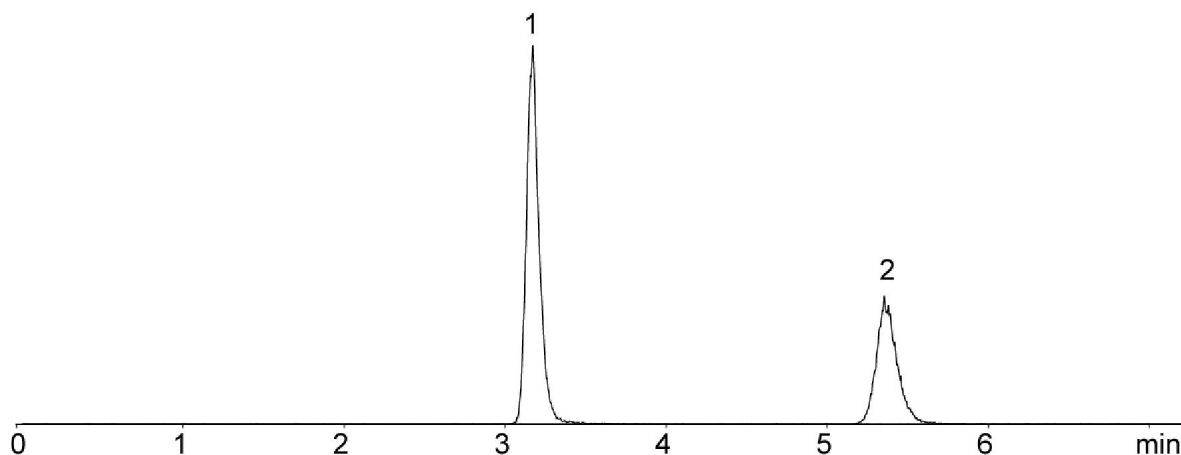

ANALYTES:

- 1 Ochratoxin B
- 2 Ochratoxin A



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.