

Tryptic Digests of three genetic variants of cytochrome c proteins on Jupiter Proteo LC column

Page 1 of 1

Column: Jupiter 4u Proteo 90A

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4396-E0

Elution Type: Gradient

Eluent A: 0.12% TFA in water

Eluent B: 0.10% TFA in MeCN


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	50	45	55
	3	55	5	95
	4	60	95	5
	5	65	95	5

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (25 °C)

Products used in this Application:

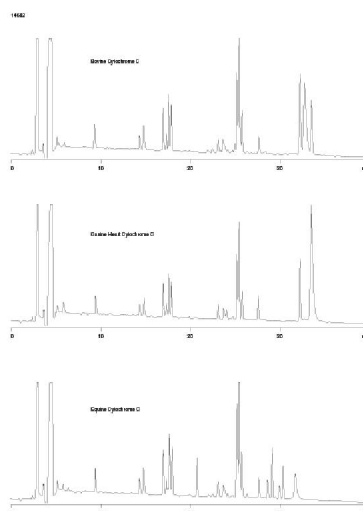

ANALYTES:

- 1 Bovine Cytochrome C
- 2 Canine Heart Cytochrome C
- 3 Equine Cytochrome C

