


## Antioxidants

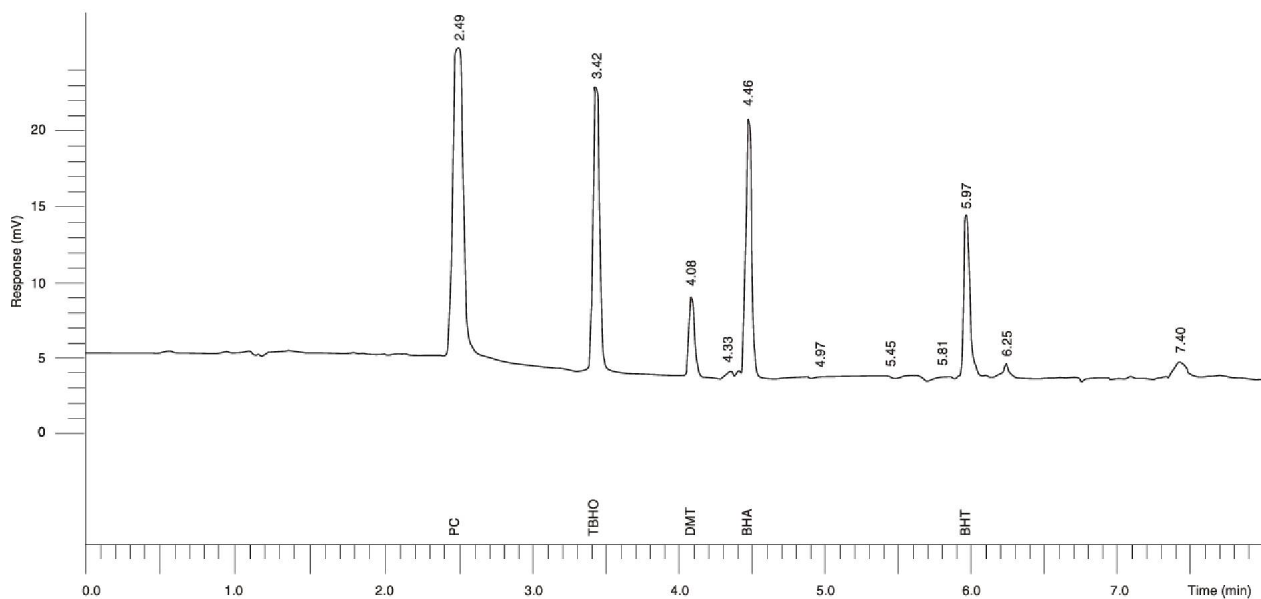

**Column:** Luna 5  $\mu$ m C18 / ODS (2)

**Dimensions:** 100 x 4.6 mm ID

**Order No:**

**Elution Type:** Gradient

**Eluent A:** Acetonitrile

**Eluent B:** 0.1% Phosphate Buffer

|                          |
|--------------------------|
| <b>Gradient Profile:</b> |
|--------------------------|

| Step No. | Step Time | Pct A | Pct B |
|----------|-----------|-------|-------|
| 1        | 0         | 30    | 70    |
| 2        | 5         | 70    | 30    |

**Flow Rate:** 1 mL/min

**Temperature:** 22 °C

**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm

**SPE:**

**SAMPLE ANALYTE**

- 1 PG
- 2 TBHQ
- 3 DMT
- 4 BHA
- 5 BHT

**For more information call your Phenomenex Representative at (310) 212-0555**

©Copyright 2000-2004 Phenomenex Inc.

**Ref. No.:** 1283/ALG/129