

Fumonisin B1 using Gemini

Column: Gemini 5u C18 110A

Dimensions: 150 x 3 mm ID

Order No: 00F-4435-Y0

Elution Type: Isocratic

Eluent A: 0.1% Formic acid in Water

Pct. A: 50

Eluent B: Acetonitrile


Pct. B: 50

Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) (ambient)


Page 1 of 1

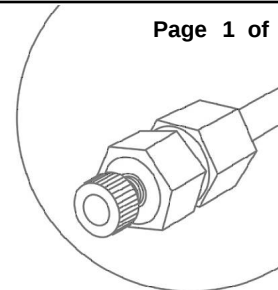

Products used in this Application:



15536

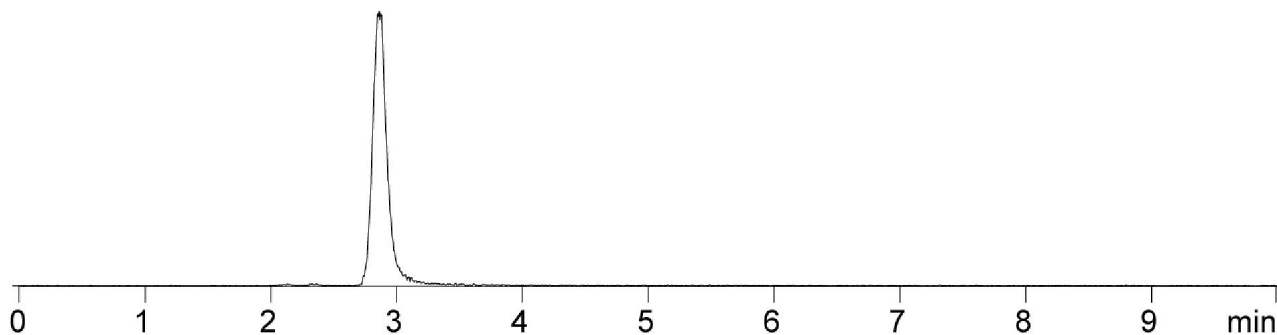

ANALYTES:

- 1 Fumonisin B1



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country,
contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.