

## Tryptic Digest of Cytochrome c on Luna SCX

Column: Luna 5u SCX 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Gradient

Eluent A: 20mM KPhos pH 2.5 + 25% Acetonitrile


Eluent B: A + 350mM KCl

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 1          | 100   | 0     |
|                   | 2        | 50         | 0     | 100   |


Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs. - Variable Wave. (ambient)

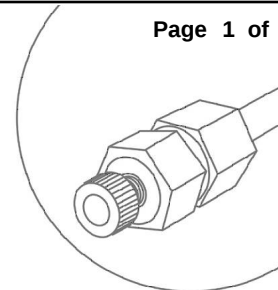

Products used in this Application:

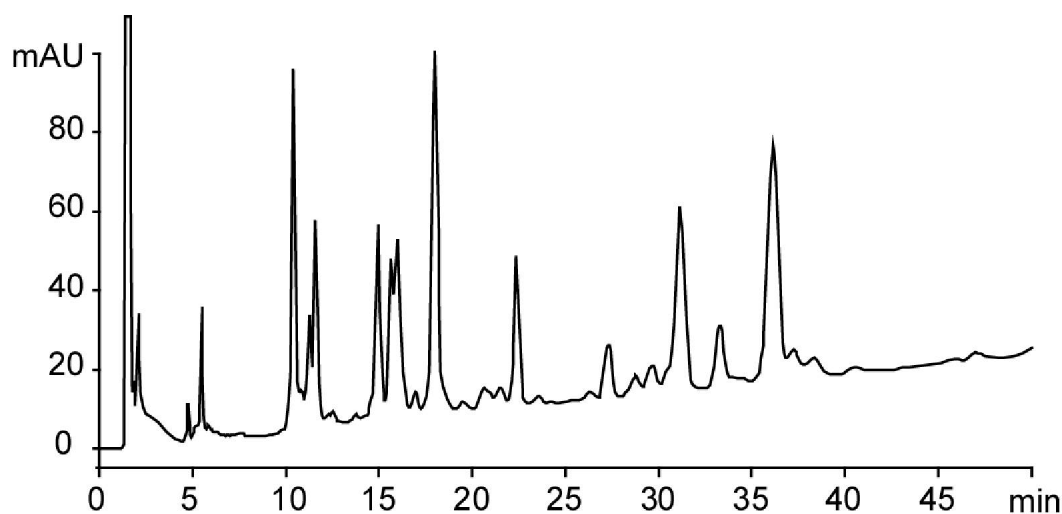**LUNA**<sup>™</sup>

## ANALYTES:

- 1 Tryptic digest of Cytochrome c



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

**USA**  
tel.: (310) 212-0555  
fax: (310) 328-7768  
email: info@phenomenex.com

**Puerto Rico**  
(800) 541-HPLC  
(310) 328-7768  
info@phenomenex.com

**Canada**  
(800) 543-3681  
(310) 328-7768  
info@phenomenex.com

**France**  
01 30 09 21 10  
01 30 09 21 11  
franceinfo@phenomenex.com

**United Kingdom**  
01625-501367  
01625-501796  
ukinfo@phenomenex.com

**Germany**  
06021-58830-0  
06021-58830-11  
anfrage@phenomenex.com

**New Zealand**  
09-4780951  
09-4780952  
info@phenomenex.co.nz

**Australia**  
02-9428-6444  
02-9428-6445  
info@phenomenex.com.au