

Drug mixture extracted from Plasma using PPT

Page 1 of 1

Column: Onyx Monolithic C18

Dimensions: 100 x 4.6 mm ID

Order No: CH0-7643

Elution Type: Gradient

Eluent A: 0.1% Formic acid in Water


Eluent B: 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	4	35	65

Flow Rate: 4 mL/min

Col. Temp.: 30 °C


Detection: Mass Spectrometer (MS) (ambient)

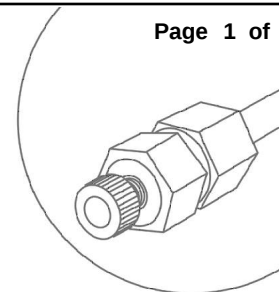

Products used in this Application:

ONYX

15587

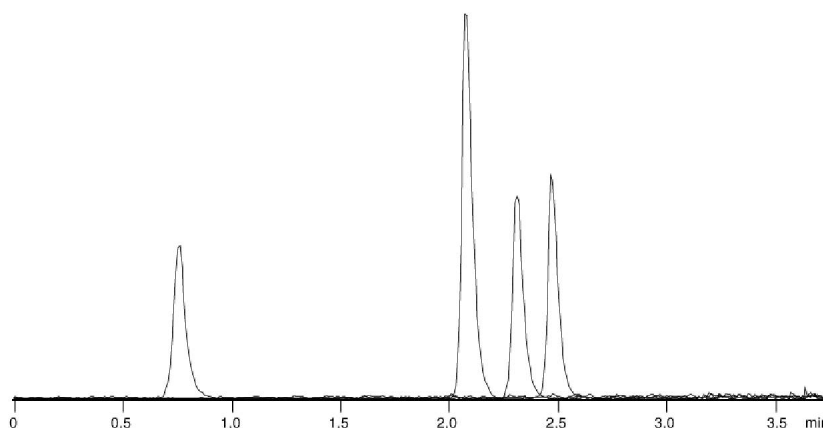

ANALYTES:

- 1 Pseudoephedrine
- 2 Propranolol
- 3 Diltiazem
- 4 Verapamil



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.