

## Carbamates

Column: Luna 3u C18(2) 100A

Dimensions: 150 x 2 mm ID

Order No: 00F-4251-B0

Elution Type: Gradient

Eluent A: Water


Eluent B: Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 90    | 10    |
|                   | 2        | 20         | 30    | 70    |
|                   | 3        | 22         | 90    | 10    |
|                   | 4        | 27         | 90    | 10    |


Flow Rate: 0.2 mL/min

Col. Temp.: 42 °C

Detection: Mass Spectrometer - API-ES + (ambient)

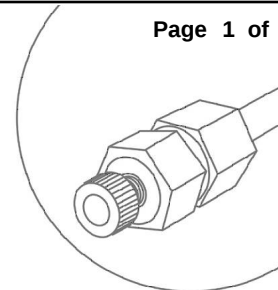

Products used in this Application:

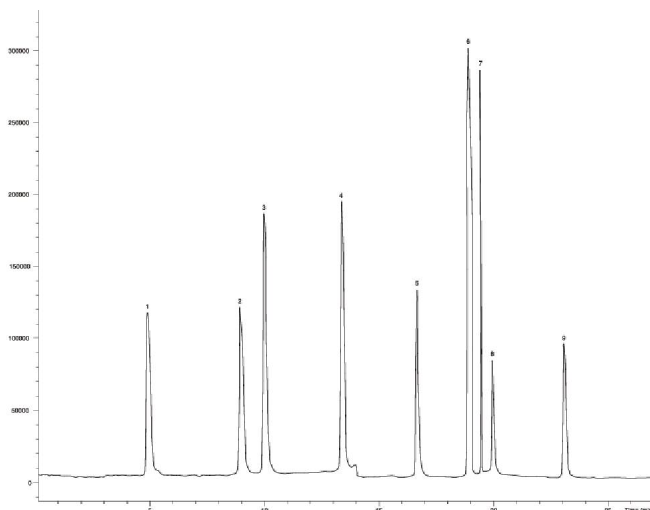**LUNA™**


## ANALYTES:

- 1 Aldicarb Sulfoxide
- 2 Oxamyl
- 3 Methomyl
- 4 3-OH-Carbofuran
- 5 Aldicarb
- 6 Propoxur
- 7 Carbofuran
- 8 Carbaryl
- 9 Methiocarb



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.