


APPLICATION NOTE - HPLC
IgG_Dog (Native) Jupiter 300 C4



AppID 14908

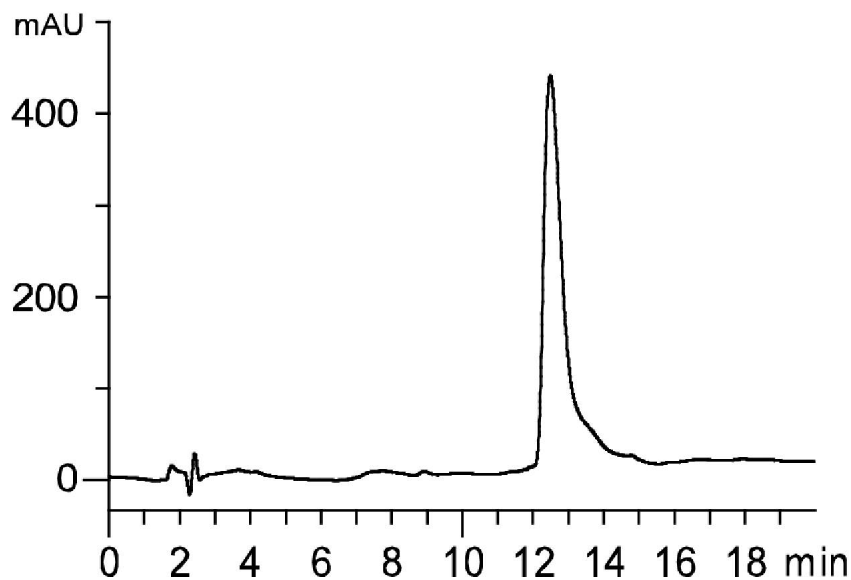

Column: Jupiter 5 μ m C4
Dimensions: 150 x 2 mm ID
Order No: 00F-4167-B0
Elution Type: Gradient
Eluent A: 0.1%TFA/95% Water/ 5% Acetonitrile
Eluent B: 0.085% TFA/75% Acetonitrile/ 20% IPA/5% Water

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	80	20
2	20	10	90

Flow Rate: 0.25 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Diode Array @ nm
SPE:

SAMPLE ANALYTES:

1 IgG Dog

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2004 Phenomenex Inc.

Ref. No.: 14908

Page 1 of 1

