


APPLICATION NOTE - HPLC
IgG Dog Reduced Jupiter C4 300



14909

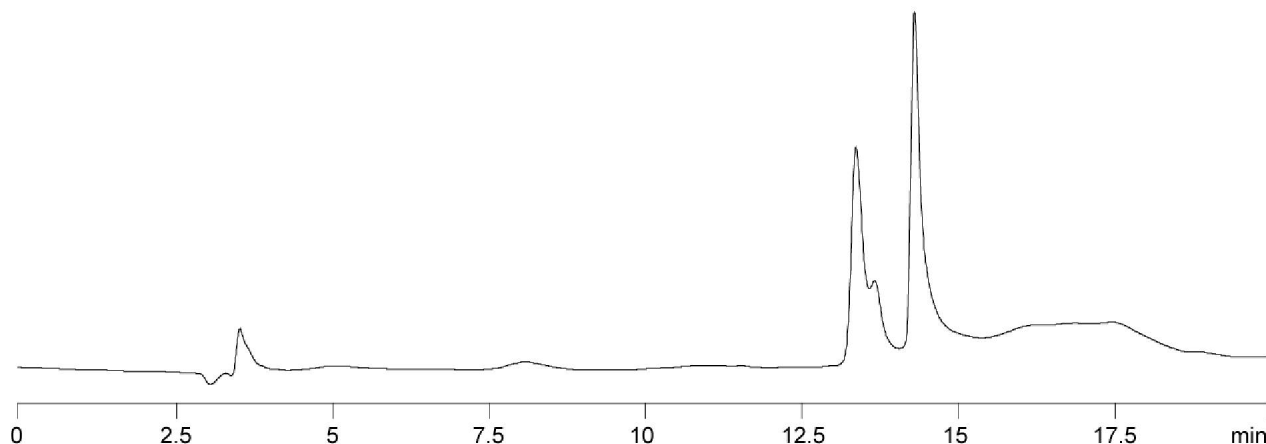

Column: Jupiter 5 μ m C4
Dimensions: 150 x 2 mm ID
Order No: 00F-4167-B0
Elution Type: Gradient
Eluent A: 0.1% TFA in Water/Acetonitrile (95:5)
Eluent B: 0.085%TFA in Acetonitrile/IPA/Water (75:20:5)

Gradient Profile:

Step No.	Step Time	Pct A	Pct B
1	0	80	20
2	20	5	95

Flow Rate: 0.25 mL/min
Temperature: 22 °C
Detection: UV-Vis Abs. - Diode Array @ 220 nm
SPE:

SAMPLE ANALYTES:

1 IgG Dog Reduced

For more information call your Phenomenex Representative at (310) 212-0555

©Copyright 2000-2004 Phenomenex Inc.

Ref. No.: 14909

Page 1 of 1

