


Alcohols

Column: Zebron ZB-5ms
Dimensions: 30meter x 0.25mm x 0.5µm
Order No.: 7HG-G010-17
Injection: Split 150:1 @ 200°C, 0.2µL
Carrier Gas: Constant Flow Helium 1.4mL/min
Oven Program: 35°C for 4 min to 100°C at 10°C/min to 250°C at 25°C/min
Detector: MSD; 20-300amu
Sample:

Gold Seal: Phenomenex single groove (AGO-7518)
Liner: Single Taper Gooseneck (AG0-4680)
Septum: Phenomenex PhenoRed-400 (AG0-4690)

- | | | |
|------------------------|--------------------------|--------------------------|
| 1. Methanol | 11. 2-Pentanol | 21. 3-Methylcyclohexanol |
| 2. Ethanol | 12. 3-Pentanol | 22. 4-Methylcyclohexanol |
| 3. Isopropanol | 13. 1-Pentanol | 23. n-Heptanol |
| 4. tert-Butanol | 14. Cyclopentanol | 24. 2-Ethylhexanol |
| 5. Allyl alcohol | 15. 3-Hexanol | 25. Benzylalcoholol |
| 6. n-Propanol | 16. 2-Ethylbutanol | 26. Phenylethanol |
| 7. 2-Butanol | 17. n-Hexanol | 27. n-Octanol |
| 8. 2-Methyl-1-propanol | 18. Cyclohexanol | 28. 2-Phenylethanol |
| 9. Methallyl alcohol | 19. 4-Heptanol | 29. Benzenepropanol |
| 10. n-Butanol | 20. 2-Methylcyclohexanol | 30. n-Decanol |

Notes: Equal mixture of pure alcohols.

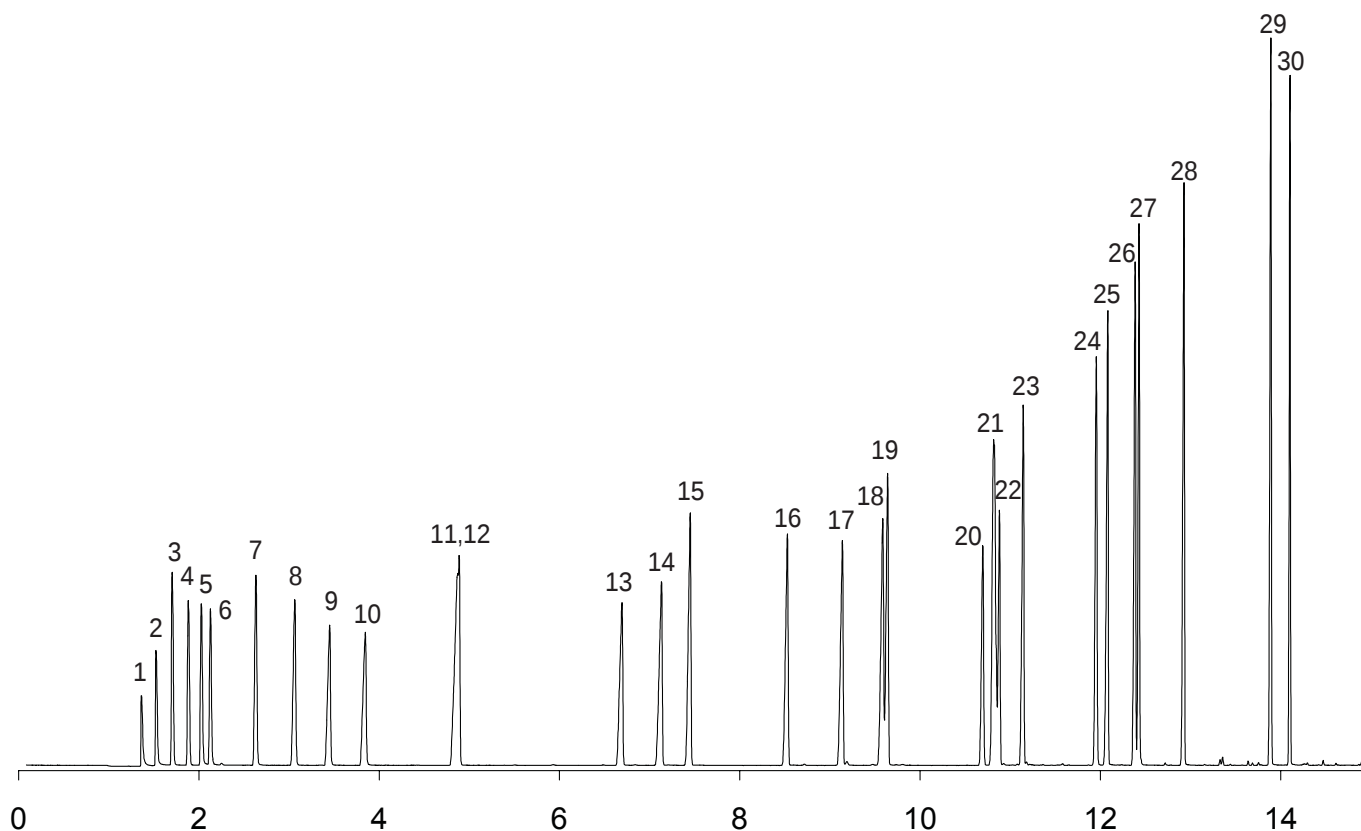

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

 **phenomenex®**
 ...breaking with tradition™

USA
 mail: 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 tel.: (310) 212-0555
 fax: (310) 328-7768
 email: info@phenomenex.com

Puerto Rico
 273 Sierra Morena,
 Suite #104
 San Juan,
 Puerto Rico 00926
 (800) 541-HPLC
 (310) 328-7768
 info@phenomenex.com

Canada
 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 (800) 543-3681
 (310) 328-7768
 info@phenomenex.com

United Kingdom
 Queens Avenue,
 Hurdfield Ind. Est.,
 Macclesfield, Cheshire
 SK10 2BN, England
 01625-501367
 01625-501796
 ukinfo@phenomenex.com

Germany
 Zeppelinstr. 5
 63741 Aschaffenburg
 Deutschland
 06021-58830-0
 06021-58830-11
 anfrage@phenomenex.com

New Zealand
 P.O. Box 31-601
 Milford
 Auckland
 New Zealand
 09-4780951
 09-4780952
 info@phenomenex.co.nz

Australia
 P.O. Box 4084
 Lane Cove, NSW 1595
 Australia
 02-9428-6444
 02-9428-6445
 info@phenomenex.com.au