


Aldehydes

Column: Zebtron ZB-WAX
Dimensions: 30m x 0.25mm x 0.25µm
Order No.: 7HG-G007-11
Injection: Split 20:1 @ 250°C, 1µL
Carrier Gas: Constant Flow Helium 0.8 mL/min
Oven Program: 40°C for 6 min to 70°C at 10°C/min to 230°C at 20°C/min for 5 min.
Detector: MS @ 280°C
Sample:

1. Propanal	11. Hexanal
2. Cyclopentane	12. Paraldehyde
3. Propenal	13. Heptanal
4. Butanal	14. Octanal
5. 2-Methyl-2-Propenal	15. contaminant
6. 2-Methylbutanal	16. Nonanal
7. Ethanol	17. Decanal
8. Pentanal	
9. 2-Methylpentanal	
10. Chloroform	

Notes: Analytes are 100-200ppm in pentane and chloroform.

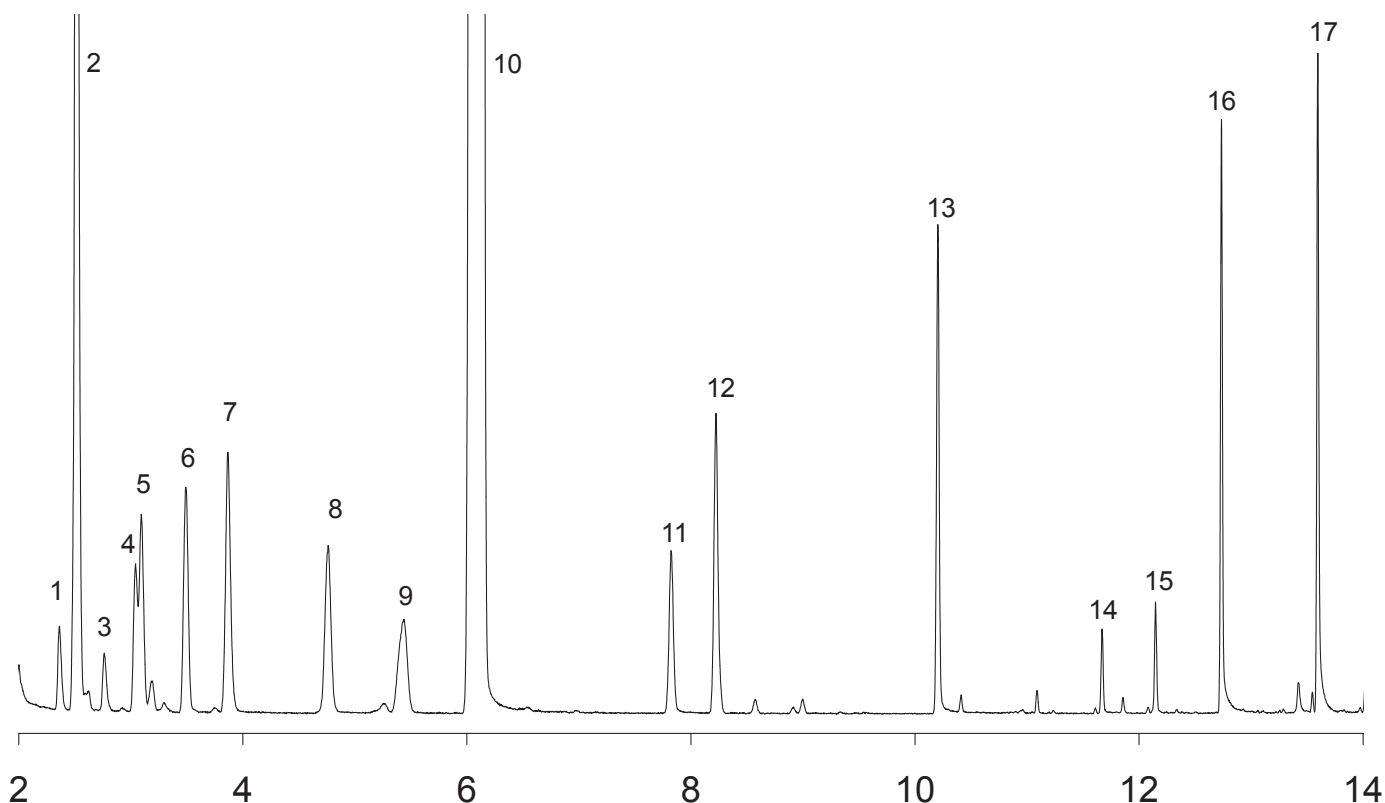

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


phenomenex®
 ...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdsfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
---	---	--	--	---	--	--