


Organic acids in wine using Rezex RHM


Column: Rezex RHM-Monosaccharide H+(8%)
Dimensions: 300 x 7.8 mm ID
Order No: 00H-0132-K0
Elution Type: Isocratic
Eluent A: 0.005N H₂SO₄
Flow Rate: 0.5 mL/min
Col. Temp.: 40 °C
Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Pct. A: 100

Page 1 of 1

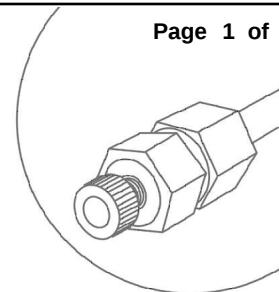

Products used in this Application:

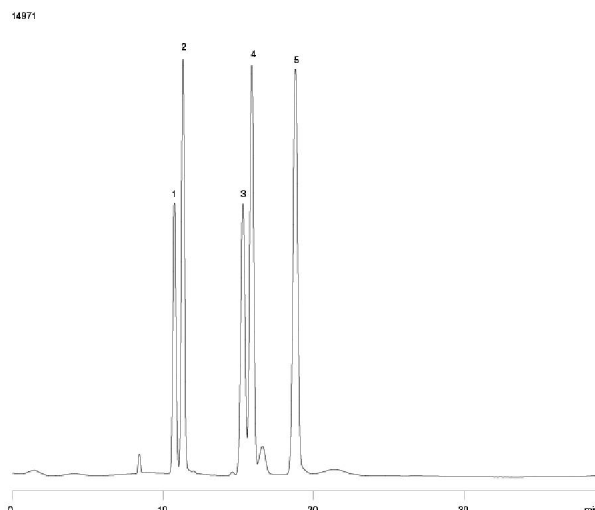REZEX[™]

ANALYTES:

- 1 Citric acid
- 2 Tartaric acid
- 3 Succinic acid
- 4 Lactic acid
- 5 Acetic acid

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.