

Sugar analysis using Rezex RCM

Column: Rezex RCM-Monosaccharide Ca+ (8%)

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0130-K0

Elution Type: Isocratic

Eluent A: Water


Pct. A: 100

Flow Rate: 0.6 mL/min

Col. Temp.: 80 °C


Detection: Evaporative Light Scattering (ambient)

Page 1 of 1

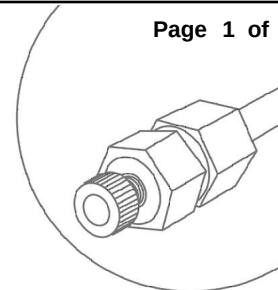

Products used in this Application:

REZEX™

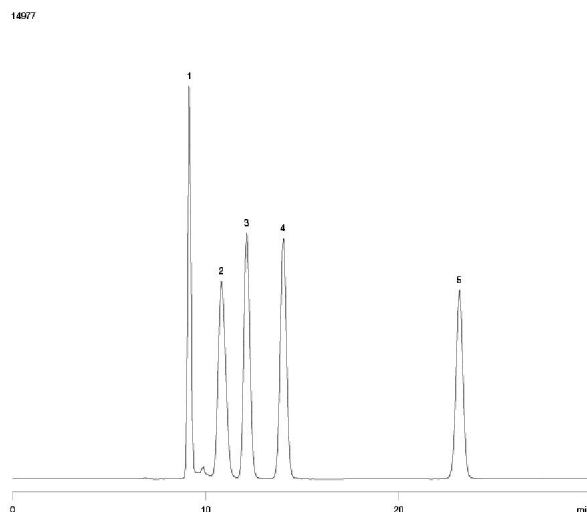

ANALYTES:

- 1 Sucrose
- 2 Glucose
- 3 Galactose
- 4 Fructose
- 5 Sorbitol



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.