


Glucosamine

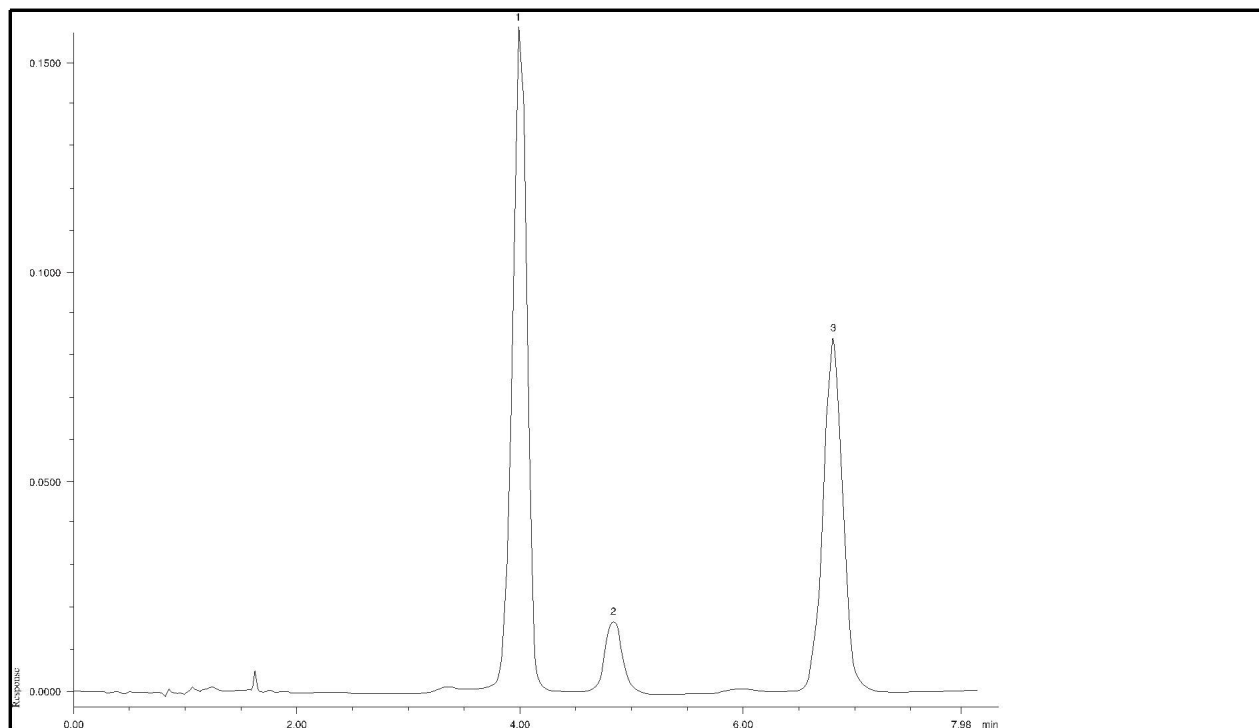

Column: Luna 5 μ m C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 10

Eluent B: Water with 0.1% H₃PO₄

Pct. B: 90

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 Glucosamine peak 1
- 2 Glucosamine peak 2
- 3 Caffeine

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 2566/ALG/176