


ASTM-2887: Simulated Distillation Markers

Column:	Zebron ZB-1																					
Dimensions:	10meter x 0.53mm x 2.65µm																					
Order No.:	7CK-G001-35																					
Injection:	Split 10:1 @ 360°C, 1µL																					
Carrier Gas:	Constant Flow Helium 8.1 mL/min																					
Oven Program:	35°C to 360°C at 15°C/min for 5 min.																					
Detector:	FID @ 360°C																					
Sample:	<table><tr><td>1. Pentane (C5)</td><td>11. Hexadecane (C16)</td></tr><tr><td>2. Hexane (C6)</td><td>12. Heptadecane (C17)</td></tr><tr><td>3. Heptane (C7)</td><td>13. Octadecane (C18)</td></tr><tr><td>4. Octane (C8)</td><td>14. Eicosane (C20)</td></tr><tr><td>5. Nonane (C9)</td><td>15. Tetracosane (C24)</td></tr><tr><td>6. Decane (C10)</td><td>16. Octacosne (C28)</td></tr><tr><td>7. Undecane (C11)</td><td>17. Dotriacontane (C32)</td></tr><tr><td>8. Dodecane (C12)</td><td>18. Hexatriacontane (C36)</td></tr><tr><td>9. Tetradecane (C14)</td><td>19. Tetracontane (C40)</td></tr><tr><td>10. Pentadecane (C15)</td><td>20. Tetratetracontane (C44)</td></tr></table>		1. Pentane (C5)	11. Hexadecane (C16)	2. Hexane (C6)	12. Heptadecane (C17)	3. Heptane (C7)	13. Octadecane (C18)	4. Octane (C8)	14. Eicosane (C20)	5. Nonane (C9)	15. Tetracosane (C24)	6. Decane (C10)	16. Octacosne (C28)	7. Undecane (C11)	17. Dotriacontane (C32)	8. Dodecane (C12)	18. Hexatriacontane (C36)	9. Tetradecane (C14)	19. Tetracontane (C40)	10. Pentadecane (C15)	20. Tetratetracontane (C44)
1. Pentane (C5)	11. Hexadecane (C16)																					
2. Hexane (C6)	12. Heptadecane (C17)																					
3. Heptane (C7)	13. Octadecane (C18)																					
4. Octane (C8)	14. Eicosane (C20)																					
5. Nonane (C9)	15. Tetracosane (C24)																					
6. Decane (C10)	16. Octacosne (C28)																					
7. Undecane (C11)	17. Dotriacontane (C32)																					
8. Dodecane (C12)	18. Hexatriacontane (C36)																					
9. Tetradecane (C14)	19. Tetracontane (C40)																					
10. Pentadecane (C15)	20. Tetratetracontane (C44)																					

Notes: Analytes are 0.1% in carbon disulfide.

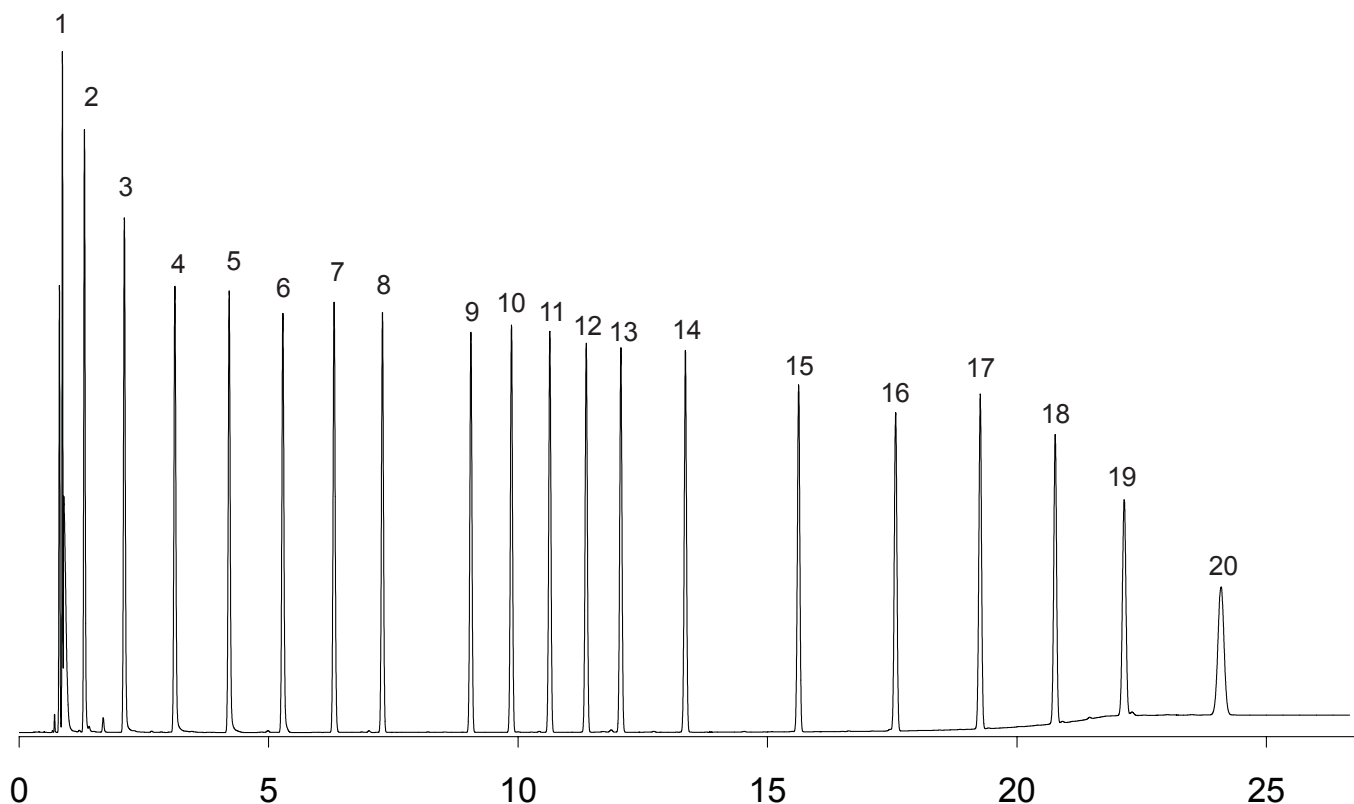

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA mail: 411 Madrid Ave. Torrance, CA 90501-1430 USA tel.: (310) 212-0555 fax: (310) 328-7768 email: info@phenomenex.com	Puerto Rico 271 Sierra Morena, Suite #104 San Juan, Puerto Rico 00926 (800) 541-HPLC (310) 328-7768 info@phenomenex.com	Canada 411 Madrid Ave. Torrance, CA 90501-1430 USA (800) 543-3681 (310) 328-7768 info@phenomenex.com	United Kingdom Queens Avenue, Hurdfield Ind. Est., Macclesfield, Cheshire SK10 2BN, England 01625-501367 01625-501796 ukinfo@phenomenex.com	Germany Zeppelinstr. 5 63741 Aschaffenburg Deutschland 06021-58830-0 06021-58830-11 anfrage@phenomenex.com	New Zealand P.O. Box 31-601 Milford Auckland New Zealand 09-4780951 09-4780952 info@phenomenex.co.nz	Australia P.O. Box 264 Pennant Hills NSW 1715 Australia 1800-553-929 1800-553-923 info@phenomenex.com.au
--	--	---	--	---	---	--