

## Onyx Si Test Chromatogram (AL0-7835)

**Column:** Onyx Monolithic Si

**Dimensions:** 100 x 4.6 mm ID

**Order No:** CH0-7648

**Elution Type:** Isocratic

**Eluent A:** n-Heptane / Dioxane

**Pct. A:** 95


**Eluent B:** Dioxane

**Pct. B:** 5

**Flow Rate:** 2 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

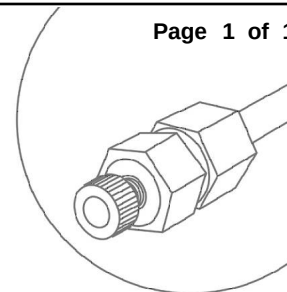


Page 1 of 1

Products used in this Application:

**ONYX**

**Analyst Note:** This standard is used to test the following Phenomenex Onyx™ Monolithic Normal Phase Columns: Onyx Si.  
NOTE: Separation conditions and results may vary for your column. Contact your Technical Representative.  
Efficiency (based on 2-Nitroanisole) = 116363 p/m (Test Mix Order No.: 7835)

15639

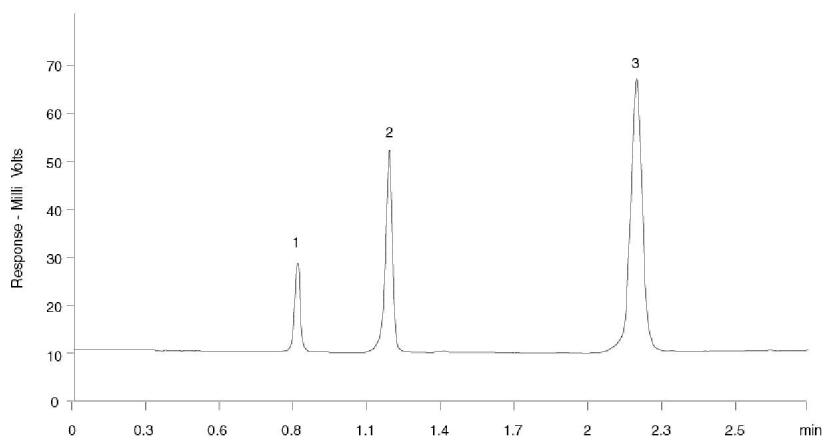

### ANALYTES:

- 1 Toluene
- 2 Nitrobenzene
- 3 2-Nitroanisole

### AUTHORIZED DISTRIBUTOR



SUNGMOON SYSTECH CORP.

#401 DaeChang Bldg., 132-9 JoongGok-Dong • KwangJin-Gu, Seoul 143-220 • Korea  
Tel: +82-2-2201-5880 • Fax: +82-2-2201-5796 • Email: info@phenomenex.co.kr