


Organic Acid screen

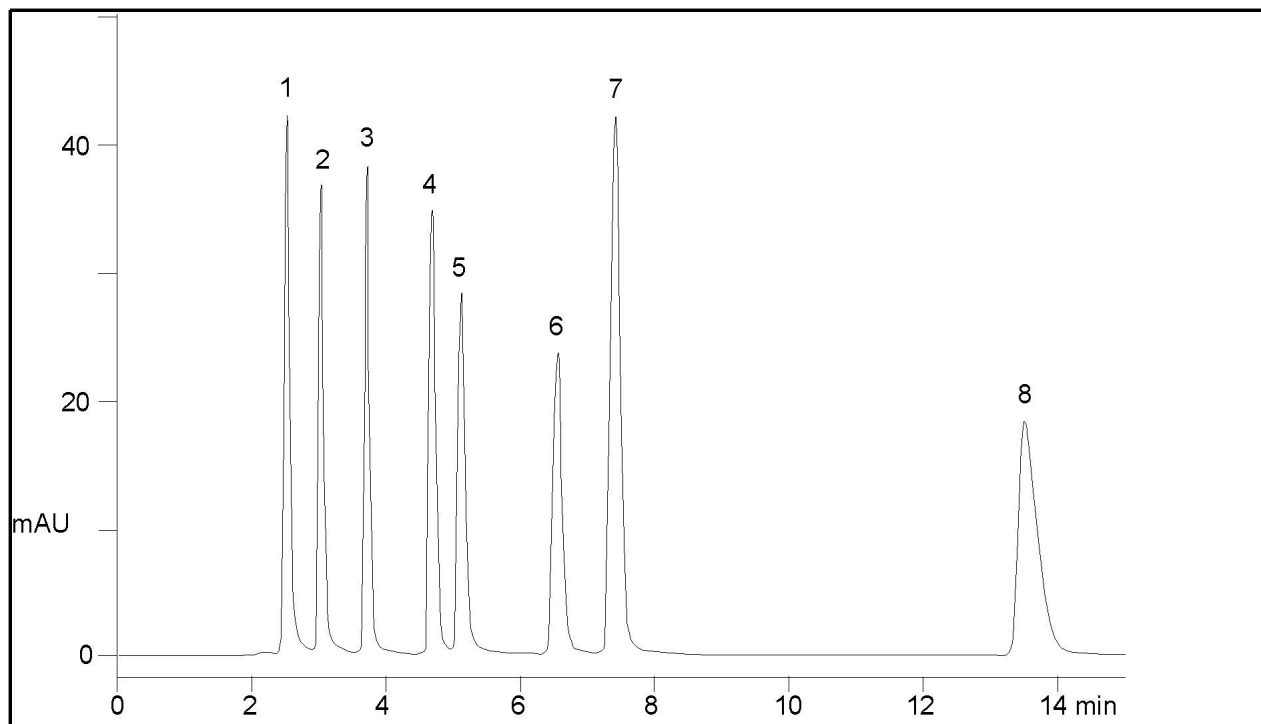

Column: Aqua 5 μ m C18 / ODS 125Å

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate

Pct. A: 97

Eluent B: MeOH

Pct. B: 3

Flow Rate: 0.7 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm

SAMPLE ANALYTE

- 1 Oxalic acid
- 2 Glycolic acid
- 3 Malic acid
- 4 Lactic acid
- 5 Acetic acid
- 6 Citric acid
- 7 Fumaric acid
- 8 Propionic acid

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 3306