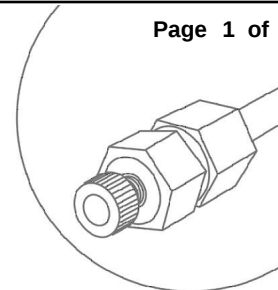
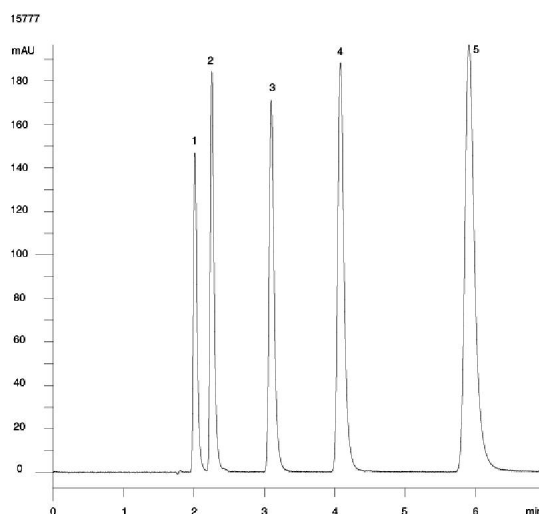


## Nucleobases - Onyx Reproducibility (Batch D)

**Column:** Onyx Monolithic C18**Dimensions:** 100 x 4.6 mm ID**Order No:** CH0-7643**Elution Type:** Isocratic**Eluent A:** Methanol**Pct. A:** 8**Eluent B:** 20mM Ammonium acetate, pH 4.5**Pct. B:** 92**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)**Comment:** Batch D

Page 1 of 1

**Products used in this Application:****ONYX****Analyst Note:** Batch D**ANALYTES:**

- 1 Cytosine
- 2 Uracil
- 3 Guanine
- 4 Thymine
- 5 Adenine

**www.phenomenex.com**

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.