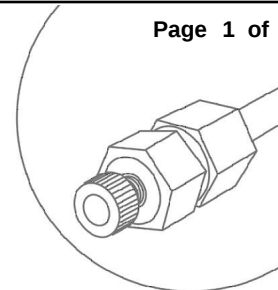
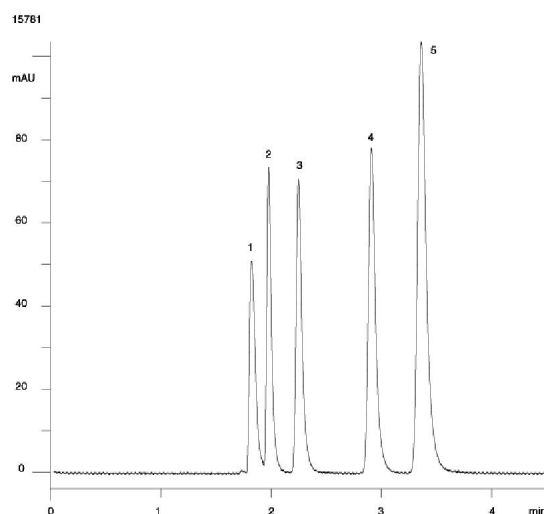


Nucleosides - Onyx Reproducibility (Batch D)

Column: Onyx Monolithic C18**Dimensions:** 100 x 4.6 mm ID**Order No:** CH0-7643**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 5**Eluent B:** 20mM Ammonium acetate, pH 4.5**Pct. B:** 95**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 260 nm (ambient)**Comment:** Batch D

Page 1 of 1

**Products used in this Application:**ONYX[™]**Analyst Note:** Batch D**ANALYTES:**

- 1 Cytidine
- 2 Uridine
- 3 Guanosine
- 4 Thymidine
- 5 Adenosine

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.