


FAME Analysis

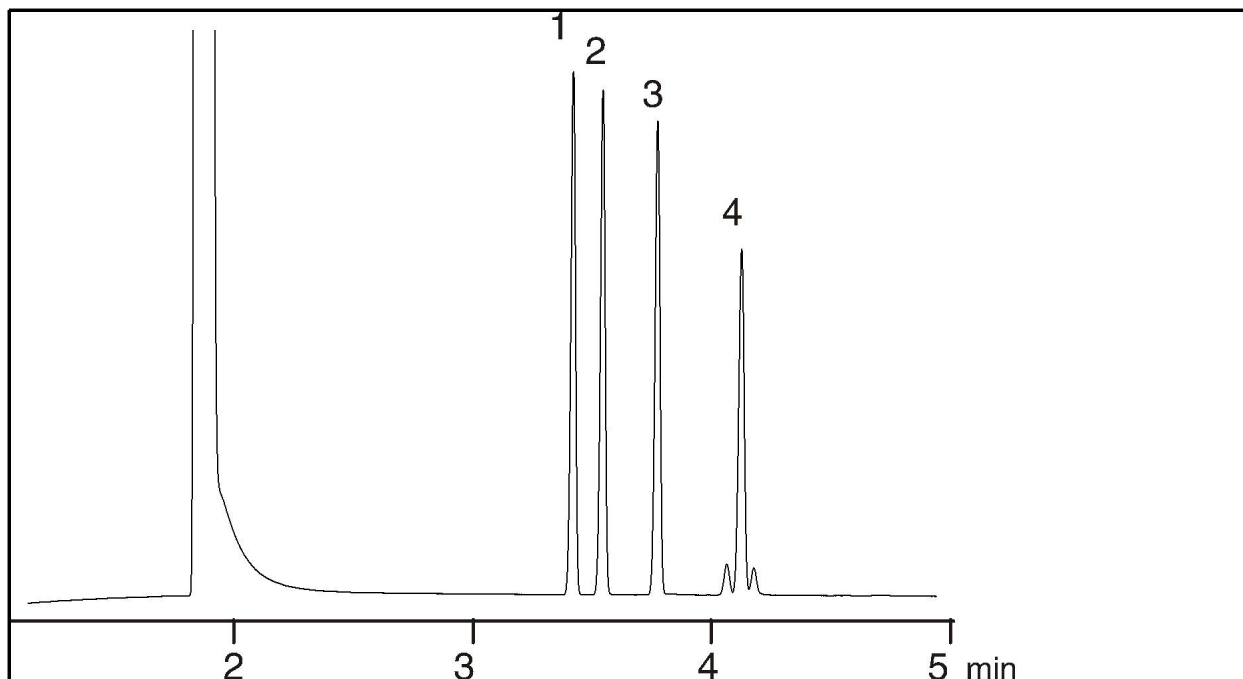

Column: Zebron ZB-WAX
Phase: Polyethylene glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 μ m
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	0	205
2	1	0	225
3	6	0	225

Carrier Gas: Helium 30 cm/sec @ 225 °C
Flow Rate: 0.88 mL/min
Injection: Split 50:1 1 μ L @ 250 °C
Detection: Flame Ionization (FID) 250 °C

SAMPLE ANALYTE

- 1 Methyl stearate
- 2 Methyl oleate
- 3 Methyl linoleate
- 4 Methyl linolenate

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. 1333/ZEB/07002 / GC_26b.wmf
No.: