


Ginkgo Leaf (Ginkgo biloba)

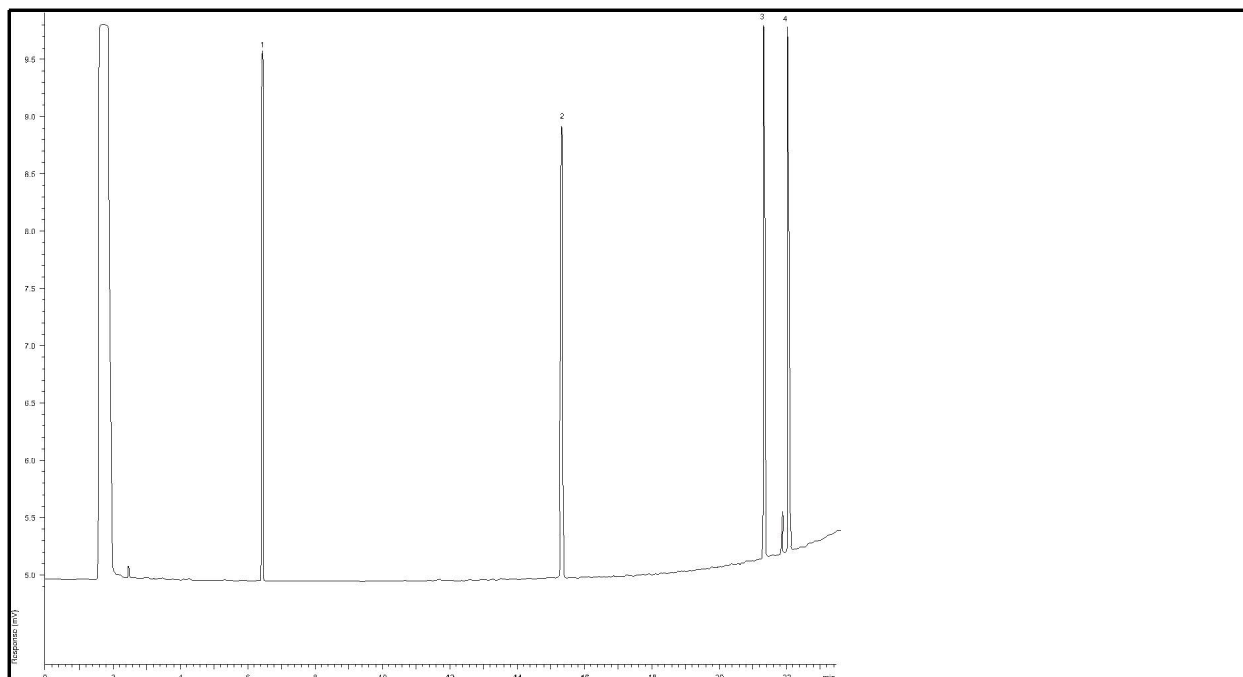

Column: Zebron ZB-1
Phase: 100% dimethyl polysiloxane
Dimensions: 30 meters x 0.32 mm x 0.25 μ m
Order No:

Gradient Profile:

Step No.	Step Time	Flow Rate	Temperature
1	0	0	120
2	24	0	340

Carrier Gas: Helium 35.23 cm/sec
Flow Rate: 1.7 mL/min
Injection: splitless 15:5 1 μ L @ 280 $^{\circ}$ C
Detection: Flame Ionization (FID) 300 $^{\circ}$ C

SAMPLE ANALYTE

- 1 Benzophenone
- 2 Bilobalide
- 3 Ginkgolide A
- 4 Ginkgolide B

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. 1122/NUT/ginkgo2.doc
No.: