



Human Serum Albumin on BioSep-SEC-S3000

Column: BioSep-SEC-S 3000**Dimensions:** 300 x 7.8 mm ID**Order No:** 00H-2146-K0**Elution Type:** Isocratic**Eluent A:** 10 mM NaH₂PO₄, pH 7.2**Pct. A:** 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 215 nm (ambient)

Page 1 of 1

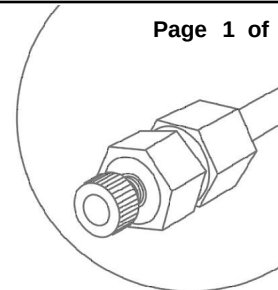

Products used in this Application:

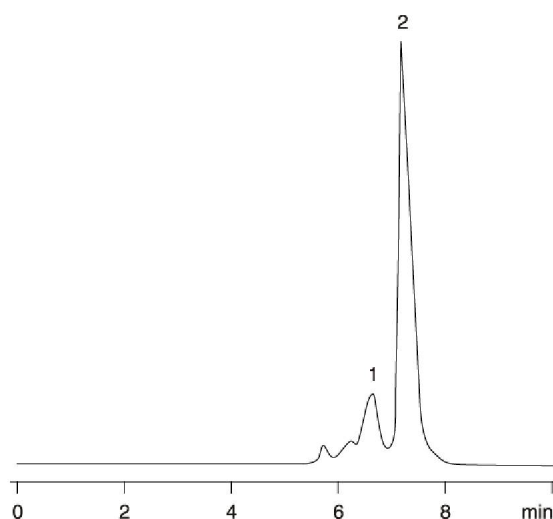BioSep[™]**ANALYTES:**

- 1 Aggregate
- 2 Albumin

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.