


Haloacetic Acids - EPA 552.2

Column: Zebron ZB-5ms
Dimensions: 30meter x 0.25mm x 0.25µm
Order No.: 7HG-G010-11
Injection: Split 20:1 @ 200°C, 1µL
Carrier Gas: Constant Flow Helium 1.9 mL/min
Oven Program: 45°C for 10 min to 145°C at 15°C/min for 5 min
Detector: ECD @ 260°C
Sample:

Gold Seal: Phenomenex single groove (AGO-7518)
Liner: Single Taper Gooseneck (AGO-4680)
Septum: Phenomenex PhenoRed-400 (AGO-4690)

1. Chloroacetic Acid
2. Bromoacetic Acid
3. Dalapon
4. Dichloroacetic Acid
5. Trichloroacetic Acid
6. Bromochloroacetic Acid
7. 1,2,3-Trichloropropane (IS)
8. Bromodichloroacetic Acid
9. Dibromoacetic Acid
10. Chlorodibromoacetic Acid
11. Tribromoacetic Acid

Notes: Internal Standard (7) was not methylated.

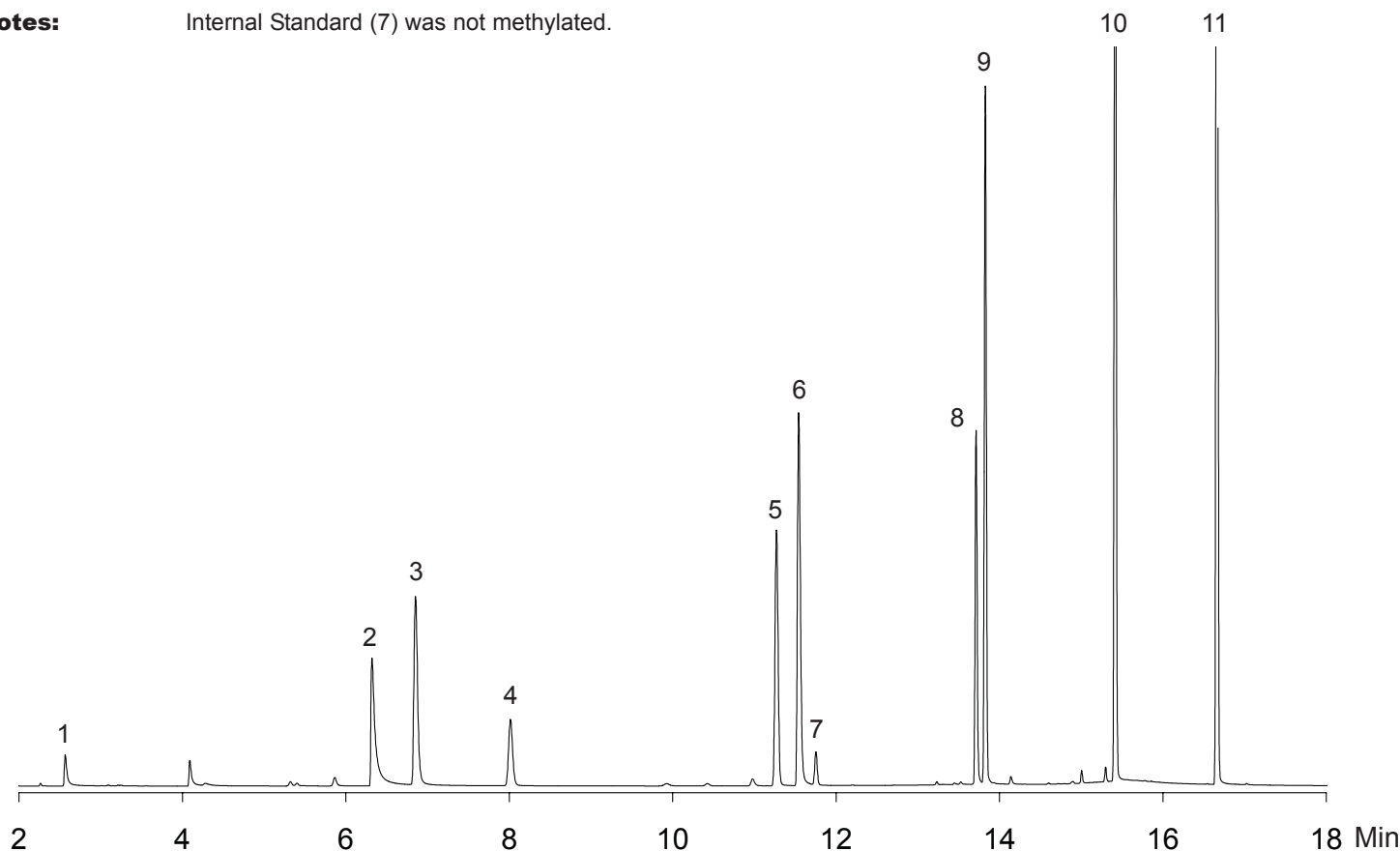

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with tradition™

USA
 mail: 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 tel.: (310) 212-0555
 fax: (310) 328-7768
 email: info@phenomenex.com

Puerto Rico
 273 Sierra Morena,
 Suite #104
 San Juan,
 Puerto Rico 00926
 (800) 541-HPLC
 (310) 328-7768
 info@phenomenex.com

Canada
 411 Madrid Avenue
 Torrance, CA
 90501-1430
 USA
 (800) 543-3681
 (310) 328-7768
 info@phenomenex.com

United Kingdom
 Queens Avenue,
 Hurdfield Ind. Est.,
 Macclesfield, Cheshire
 SK10 2BN, England
 01625-501367
 01625-501796
 ukinfo@phenomenex.com

Germany
 Zeppelinstr. 5
 63741 Aschaffenburg
 Deutschland
 06021-58830-0
 06021-58830-11
 anfrage@phenomenex.com

New Zealand
 P.O. Box 31-601
 Milford
 Auckland
 New Zealand
 09-4780951
 09-4780952
 info@phenomenex.co.nz

Australia
 P.O. Box 4084
 Lane Cove, NSW 1595
 Australia
 02-9428-6444
 02-9428-6445
 info@phenomenex.com.au