


Ketones

Column: Zebron ZB-WAX
Dimensions: 60m x 0.53mm x 1.0µm
Order No.: 7KK-G007-22
Injection: Split 10:1 @ 250°C, 1µL
Carrier Gas: Constant Flow Helium 2.9 mL/min
Oven Program: 50°C to 80°C at 10°C/min to 230°C at 25°C/min for 8 min.
Detector: FID @ 275°C
Sample:

1. Solvent Impurity	11. 5-Methyl-2-Hexanone
2. 2-Butanone	12. 3-Heptanone
3. 3-Methyl-2-Butanone	13. 2-Heptanone
4. 2-Pentanone	14. 2-Octanone
5. 3-Pentanone	15. 2-Nonanone
6. 4-Methyl-2-Pentanone	16. 2-Decanone
7. 3-Methyl-2-Pentanone	17. 2-Undecanone
8. 3-Hexanone	18. 2-Dodecanone
9. 2-Methyl-3-Hexanone	19. 2-Tridecanone
10. 2-Hexanone	

Notes: All compounds are at 100 ppm.

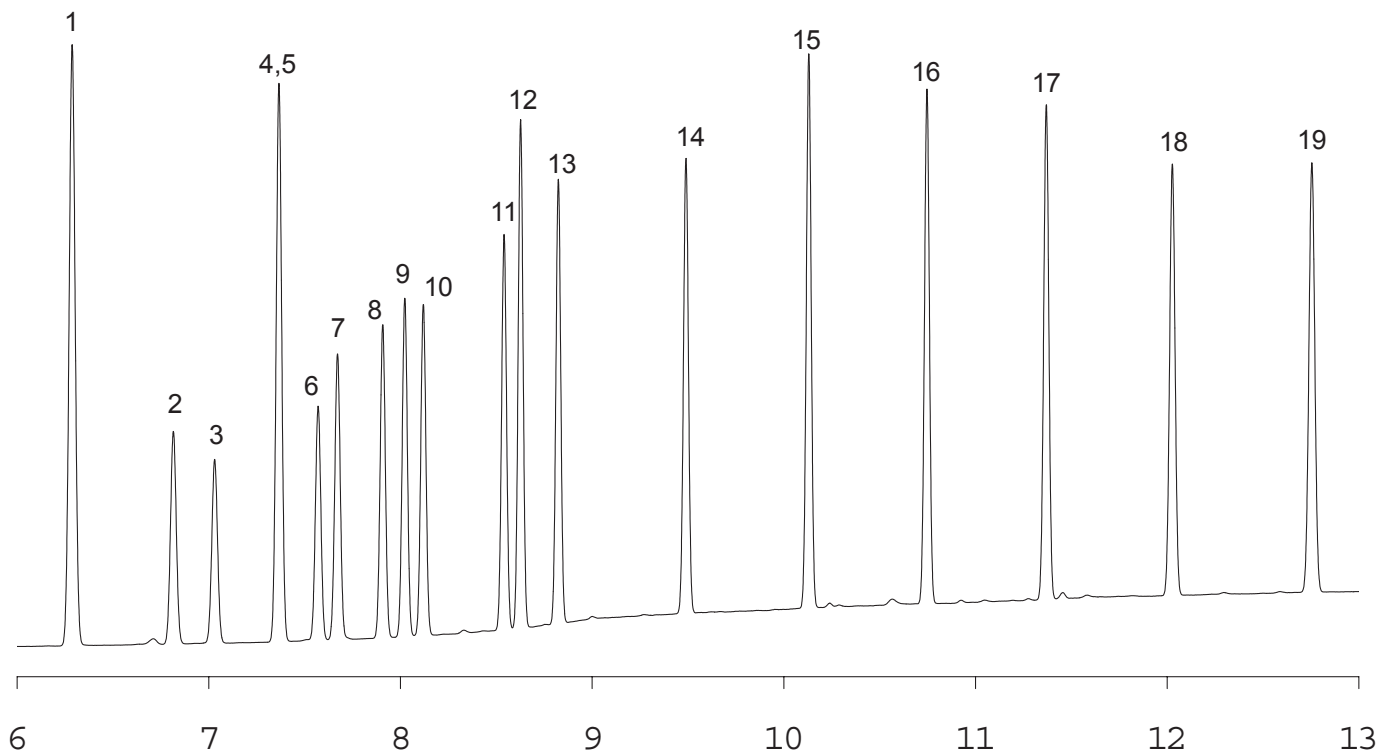

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

phenomenex®
...breaking with traditionSM

USA
 mail: 411 Madrid Ave.
 Torrance, CA
 90501-1430
 USA
 tel.: (310) 212-0555
 fax: (310) 328-7768
 email: info@phenomenex.com

Puerto Rico
 271 Sierra Morena,
 Suite #104
 San Juan,
 Puerto Rico 00926
 (800) 541-HPLC
 (310) 328-7768
info@phenomenex.com

Canada
 411 Madrid Ave.
 Torrance, CA
 90501-1430
 USA
 (800) 543-3681
 (310) 328-7768
info@phenomenex.com

United Kingdom
 Queens Avenue,
 Hurdsfield Ind. Est.,
 Macclesfield, Cheshire
 SK10 2BN, England
 01625-501367
 01625-501796
ukinfo@phenomenex.com

Germany
 Zeppelinstr. 5
 63741 Aschaffenburg
 Deutschland
 06021-58830-0
 06021-58830-11
anfrage@phenomenex.com

New Zealand
 P.O. Box 31-601
 Milford
 Auckland
 New Zealand
 09-4780951
 09-4780952
info@phenomenex.co.nz

Australia
 P.O. Box 264
 Pennant Hills
 NSW 1715
 Australia
 1800-553-929
 1800-553-923
info@phenomenex.com.au