


**USP 467 - Method I**

App ID 14797

**Column:** Zebtron ZB-5  
**Dimensions:** 30meter x 0.53mm x 5.0µm with 5m Z-guard  
**Order No.:** 7HK-G002-39; 7AK-G000-00-GZ0  
**Injection:** Split 2:1 @ 70°C, 1µL  
**Carrier Gas:** Constant Flow Helium 5.0 mL/min  
**Oven Program:** 35°C for 5 min to 175°C at 8°C/min to 260°C at 35°C/min for 16  
**Detector:** FID @ 260°C  
**Sample:**

1. Methylene Chloride
2. Chloroform
3. Benzene
4. Trichloroethylene
5. 1,4-Dioxane

**Notes:** Methylene Chloride-10ug/mL; Chloroform-1ug/mL; Benzene, 1,4- Dioxane, and Trichloroethylene at 2 ug/mL in HPLC grade water

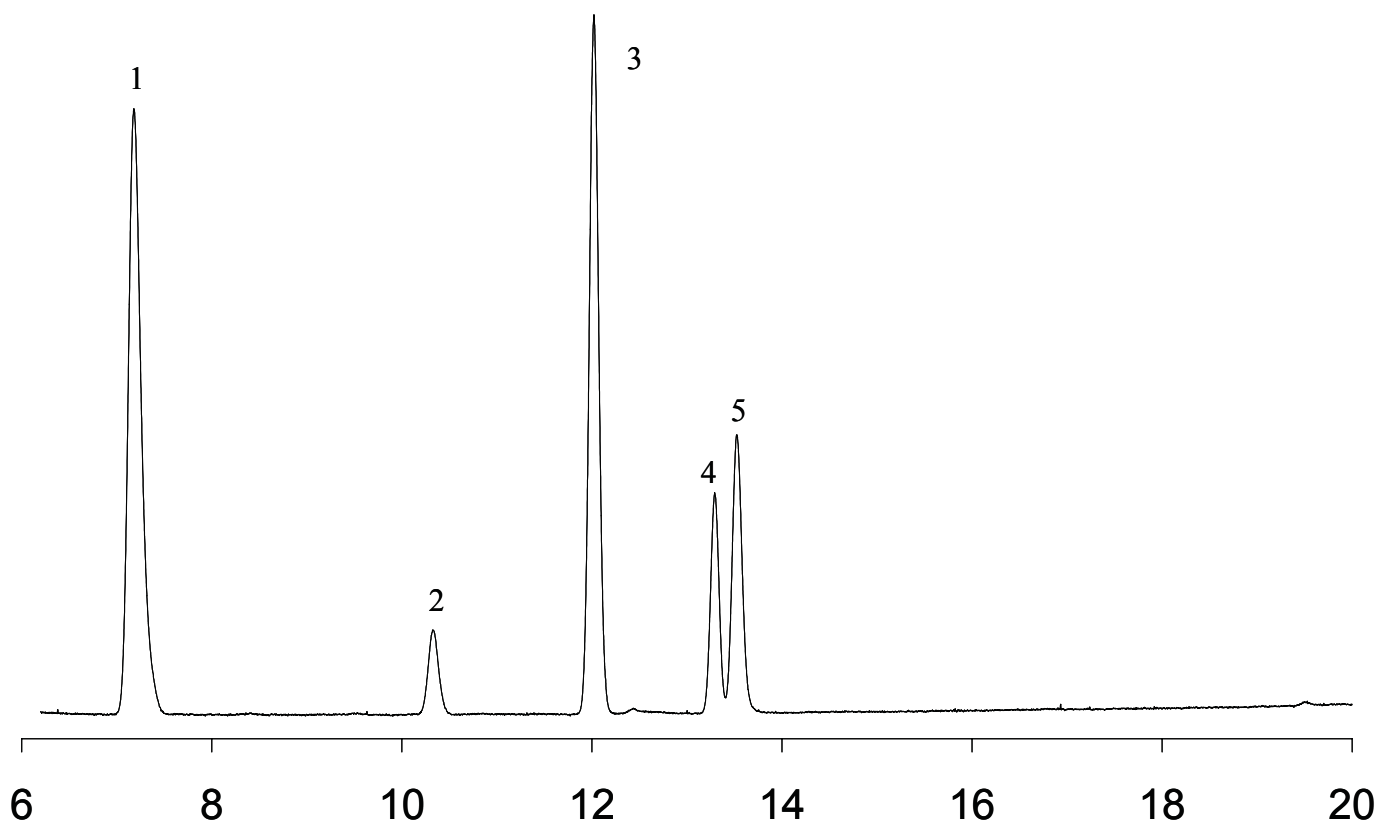

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.


**phenomenex®**  
*...breaking with tradition™*

| USA                                                          | Puerto Rico                                                        | Canada                                                 | United Kingdom                                                                        | Germany                                              | New Zealand                                           | Australia                                         |
|--------------------------------------------------------------|--------------------------------------------------------------------|--------------------------------------------------------|---------------------------------------------------------------------------------------|------------------------------------------------------|-------------------------------------------------------|---------------------------------------------------|
| mail: 411 Madrid Avenue<br>Torrance, CA<br>90501-1430<br>USA | 273 Sierra Morena,<br>Suite #104<br>San Juan,<br>Puerto Rico 00926 | 411 Madrid Avenue<br>Torrance, CA<br>90501-1430<br>USA | Queens Avenue,<br>Hurdfield Ind. Est.,<br>Macclesfield, Cheshire<br>SK10 2BN, England | Zeppelinstr. 5<br>63741 Aschaffenburg<br>Deutschland | P.O. Box 31-601<br>Milford<br>Auckland<br>New Zealand | P.O. Box 4084<br>Lane Cove, NSW 1595<br>Australia |
| tel.: (310) 212-0555                                         | (800) 541-HPLC                                                     | (800) 543-3681                                         | 01625-501367                                                                          | 06021-58830-0                                        | 09-4780951                                            | 02-9428-6444                                      |
| fax: (310) 328-7768                                          | (310) 328-7768                                                     | (310) 328-7768                                         | 01625-501796                                                                          | 06021-58830-11                                       | 09-4780952                                            | 02-9428-6445                                      |
| email: info@phenomenex.com                                   | info@phenomenex.com                                                | info@phenomenex.com                                    | ukinfo@phenomenex.com                                                                 | anfrage@phenomenex.com                               | info@phenomenex.co.nz                                 | info@phenomenex.com.au                            |