

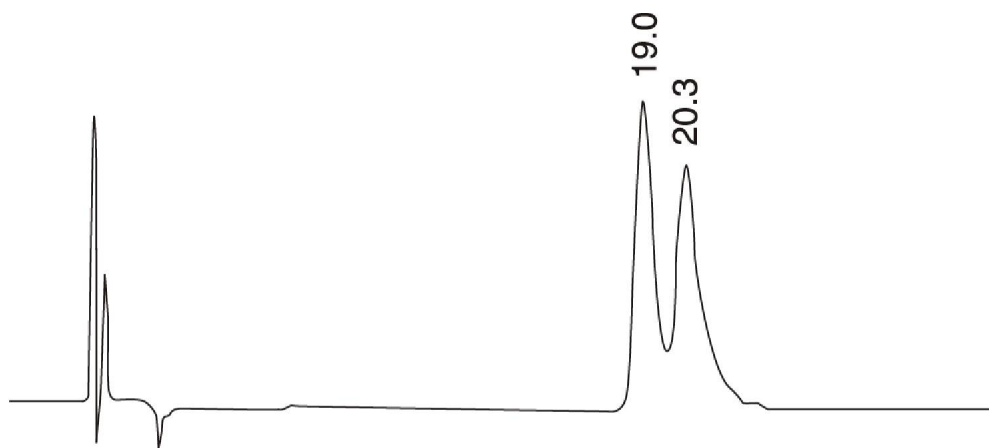
Homochlorcyclizine

Page 1 of 1

Column: Chirex 3022 (S)-indoline-2-carboxylic acid and (R)-1-(alpha-naphthyl)ethylamine
Dimensions: 250 x 4 mm ID
Order No: 00G-3022-D0
Elution Type: Isocratic
Eluent A: Hexane
Eluent B: 1,2-Dichloroethane
Eluent C: Methanol
Eluent D: Trifluoroacetic acid
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

Pct. A: 72.3**Pct. B:** 24.9**Pct. C:** 2.5**Pct. D:** 0.2

Products used in this Application:

CHIREX[®] $\alpha = 1.08$ **ANALYTES:**

- 1 Homochlorcyclizine
- 2 Homochlorcyclizine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.