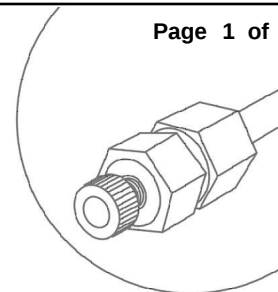
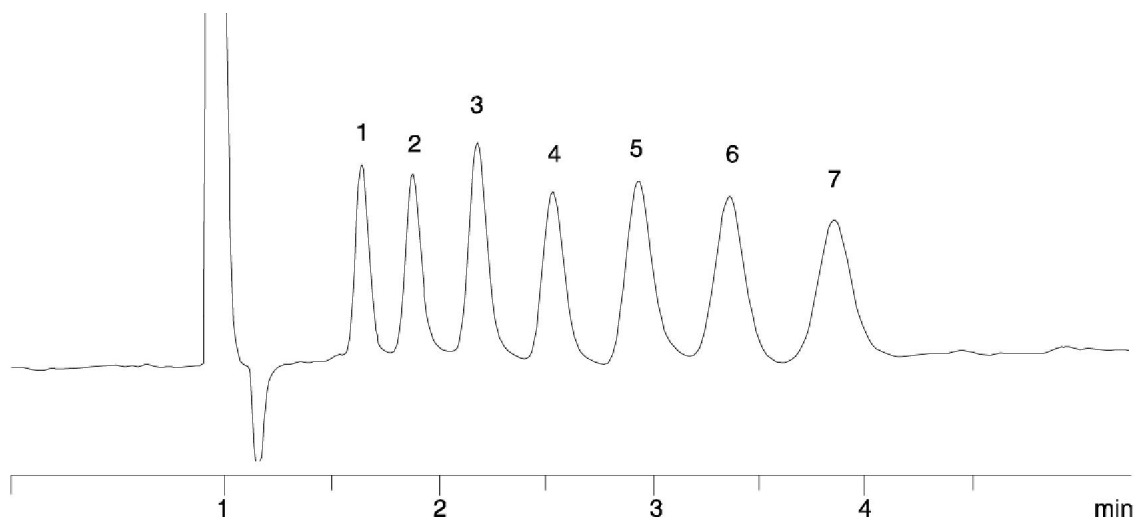


Complex Sugars on Luna Amino

Column: Luna 5u NH2 100A**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4378-E0**Elution Type:** Isocratic**Eluent A:** Acetonitrile**Pct. A:** 65**Eluent B:** Water**Pct. B:** 35**Flow Rate:** 3 mL/min**Col. Temp.:** 40 °C**Detection:** Refractive Index (RI) (40 °C)

Page 1 of 1

**Products used in this Application:****LUNA[®]****ANALYTES:**

- 1 Glucose
- 2 Maltose
- 3 Maltotriose
- 4 Maltotetraose
- 5 Maltopentaose
- 6 Maltohexaose
- 7 Maltoheptaose

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.