


Vitamins

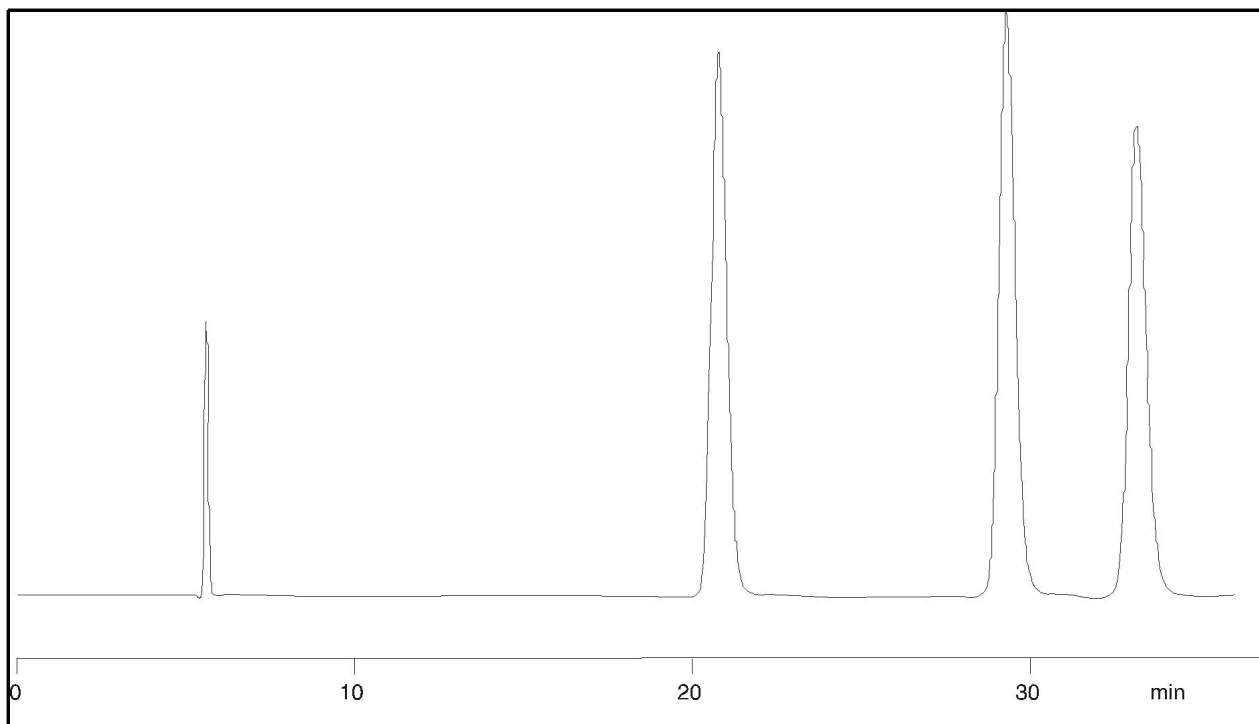

Column: Develosil 5 μ m C30 RP-Aqueous (C30-UG-5)

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 0.05M K₂HPO₄ with H₃PO₄, pH 6.9

Pct. A: 100

Flow Rate: 1 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

- 1 Orotic acid
- 2 Pyridoxal HCl
- 3 Pyridoxine HCl
- 4 Nicotinamide

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 13575/A-15