


USP Estradiol

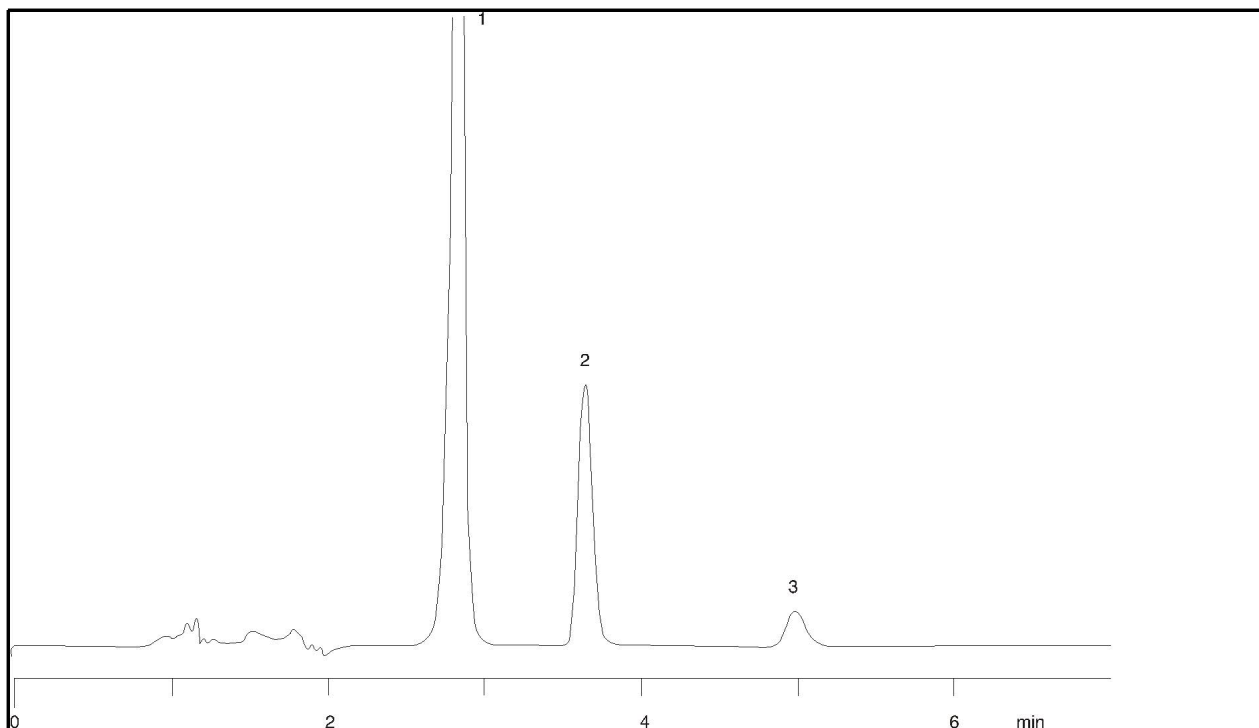

Column: Luna 5 μ m C18 / ODS (2)

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Acetonitrile

Pct. A: 55

Eluent B: Water

Pct. B: 45

Flow Rate: 1 mL/min

Temperature: 30 °C

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm

SAMPLE ANALYTE

1 *Ethyl paraben*

2 *Estrone*

3 *Estradiol*

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 10893/BRO/1