


Formic Acid as an Impurity in Acetic Acid

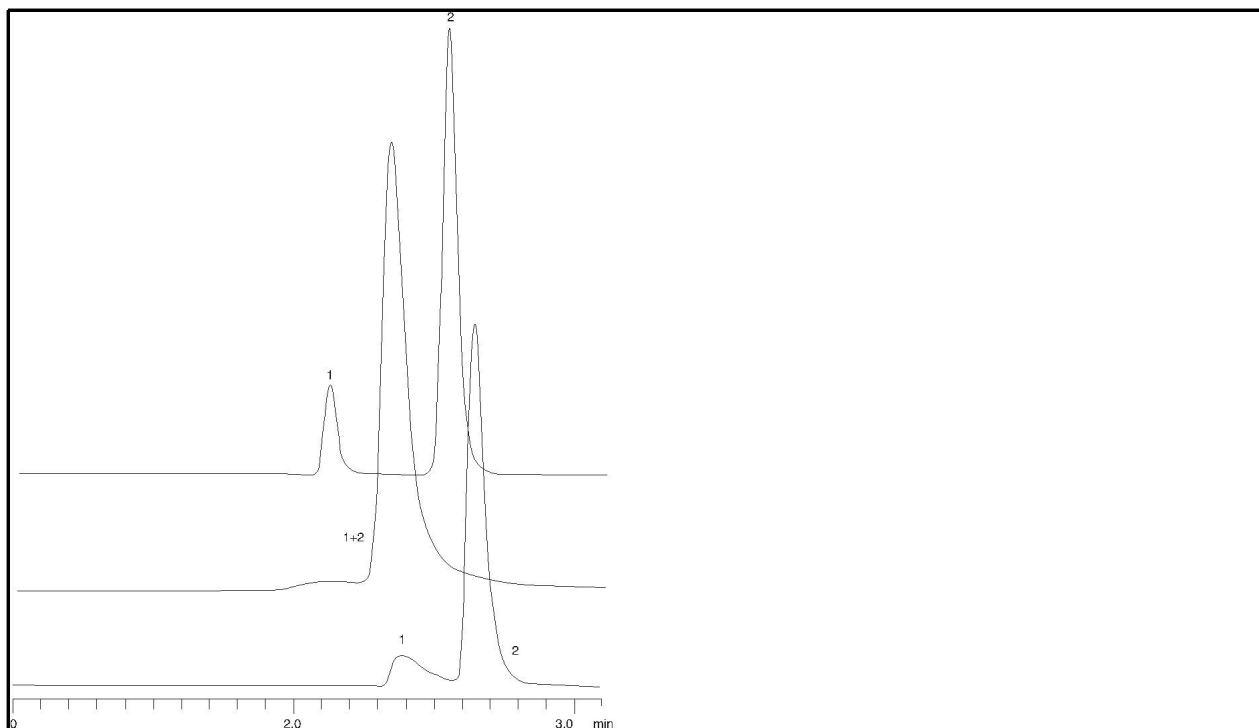

Column: Synergi 4 μ m C6H5 / Phenyl Polar-RP

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20 mM Potassium phosphate pH 2.5

Pct. A: 97

Eluent B: Methanol

Pct. B: 3

Flow Rate: 1 mL/min

Temperature: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 220 nm

SAMPLE ANALYTE

1 Formic acid

2 Acetic acid

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 10904/BRO/8