


Carbohydrates

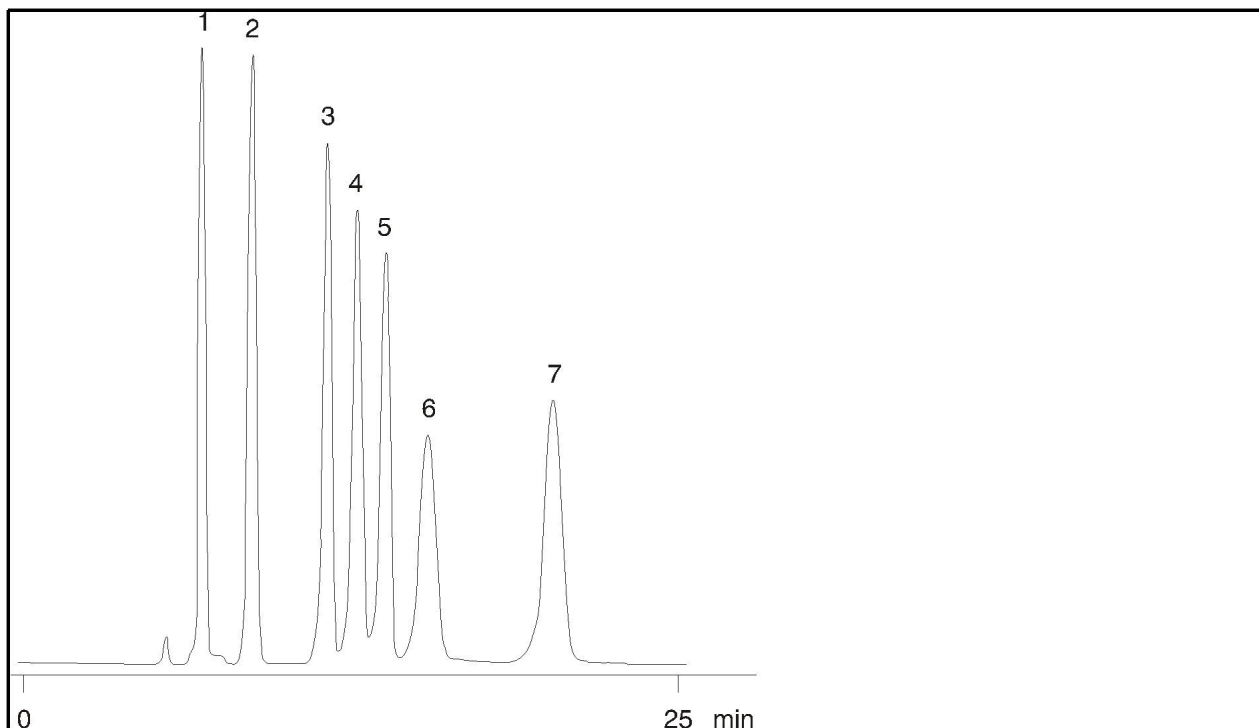

Column: Rezex 8 μ m Sulfonic Acid (SCX-SO₃H), strong cation exchange RNM-Carbohydrate

Dimensions: 300 x 7.8 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Water

Pct. A: 100

Flow Rate: 0.6 mL/min

Temperature: 75 °C

Detection: Refractive Index @ nm

SAMPLE ANALYTE

- 1 *Stachyose*
- 2 *Cellobiose*
- 3 *Glucose*
- 4 *Fructose*
- 5 *Arabinose*
- 6 *Ribose*
- 7 *Salicin (1% solution of each)*

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 5510/CAT99/236