


Ginger Root (*Zingiber officinale*)

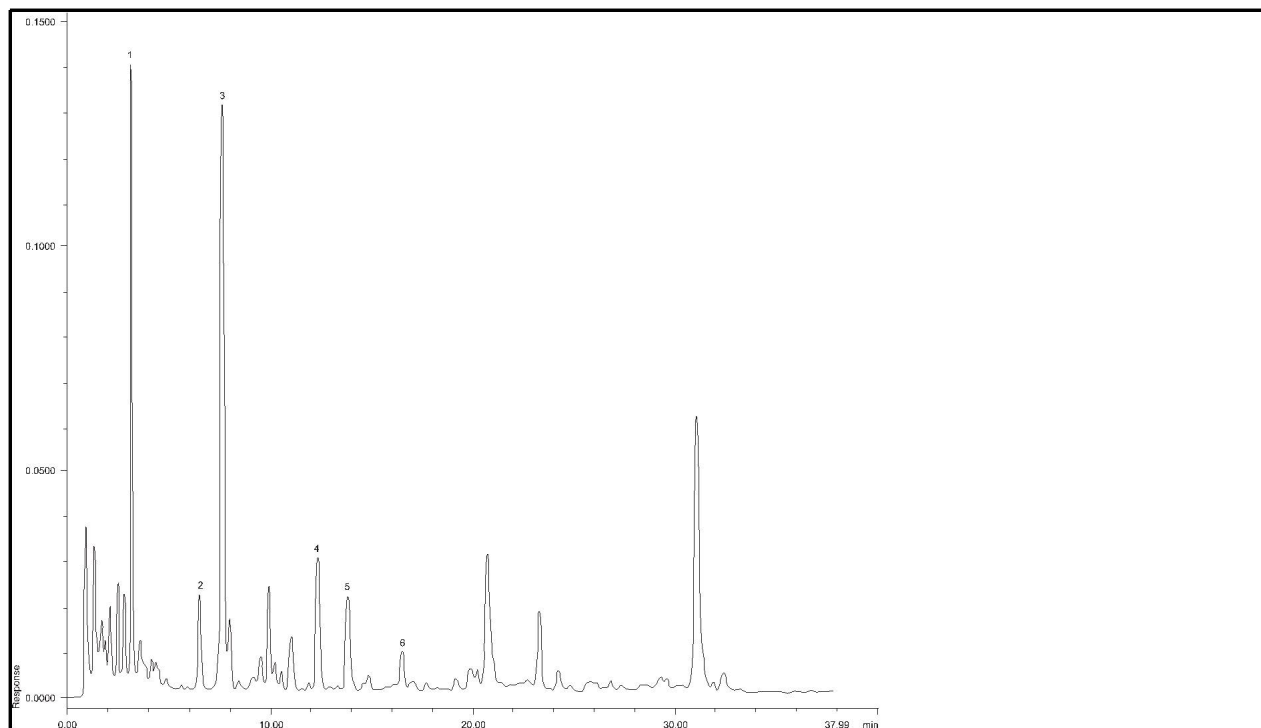

Column: Columbus 5 μ m C18 / ODS

Dimensions: 150 x 4.6 mm ID

Order No:

Elution Type: Gradient

Eluent A: Acetonitrile/water/tetrahydrofuran (50:45:5)

Eluent B: Acetonitrile/water/tetrahydrofuran (80:15:5)

Gradient Profile:		Step No.	Step Time	Pct A	Pct B
		1	0	100	0
		2	4	100	0
		3	28	0	100
		4	38	0	100

Flow Rate: 1.5 mL/min

Temperature: 25 °C

Detection: UV-Vis Abs. - Variable Wave. @ 282 nm

SPE:

SAMPLE ANALYTE

- 1 6-Gingerol
- 2 8-Gingerol
- 3 6-Shogaol
- 4 6-Gingerdiol
- 5 10-Gingerol
- 6 8-Shogaol

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1121/NUT/ginger1.doc