


Diuretics

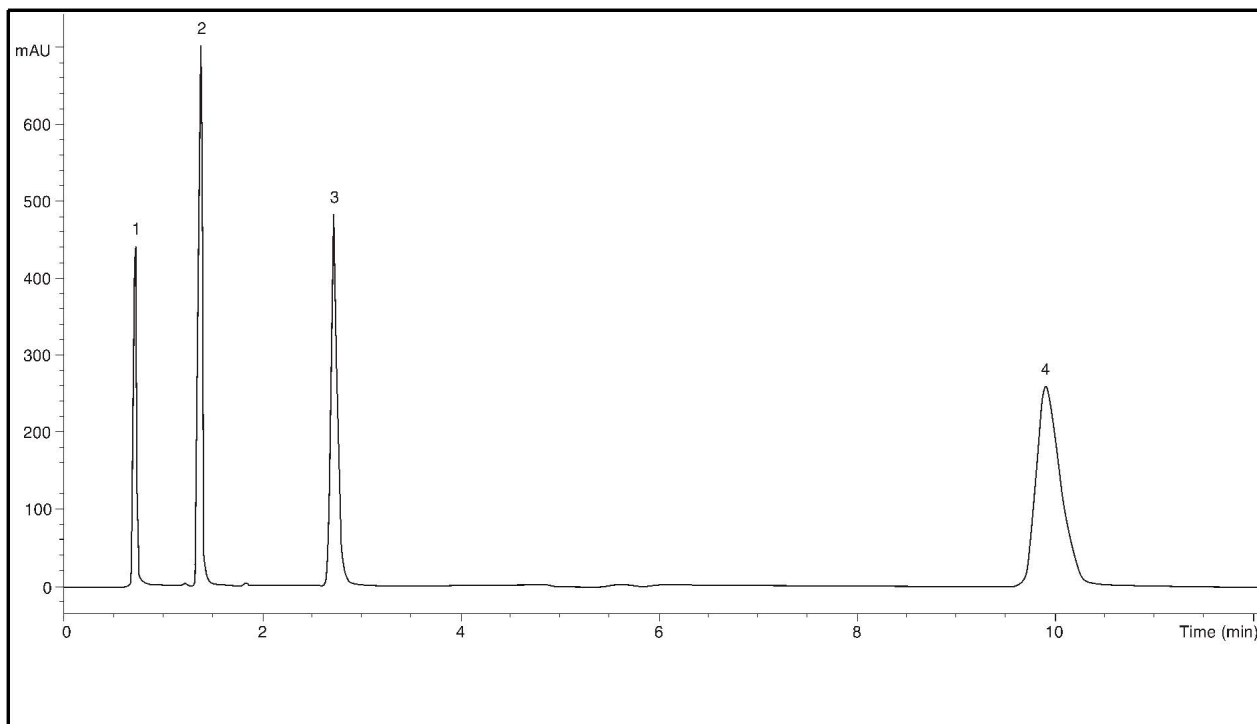

Column: Luna 3 μ m C18 / ODS (2)

Dimensions: 75 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: 20 mM Na₂HPO₄, pH 5.0

Pct. A: 74

Eluent B: Acetonitrile

Pct. B: 26

Flow Rate: 1 mL/min

Temperature: 35 °C

Detection: UV-Vis Abs. - Variable Wave. @ 270 nm

SAMPLE ANALYTE

- 1** Amiloride
- 2** Hydrochlorothiazide
- 3** Furosemide
- 4** Bumetanide

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1197/ALG/46