


Sun Screens

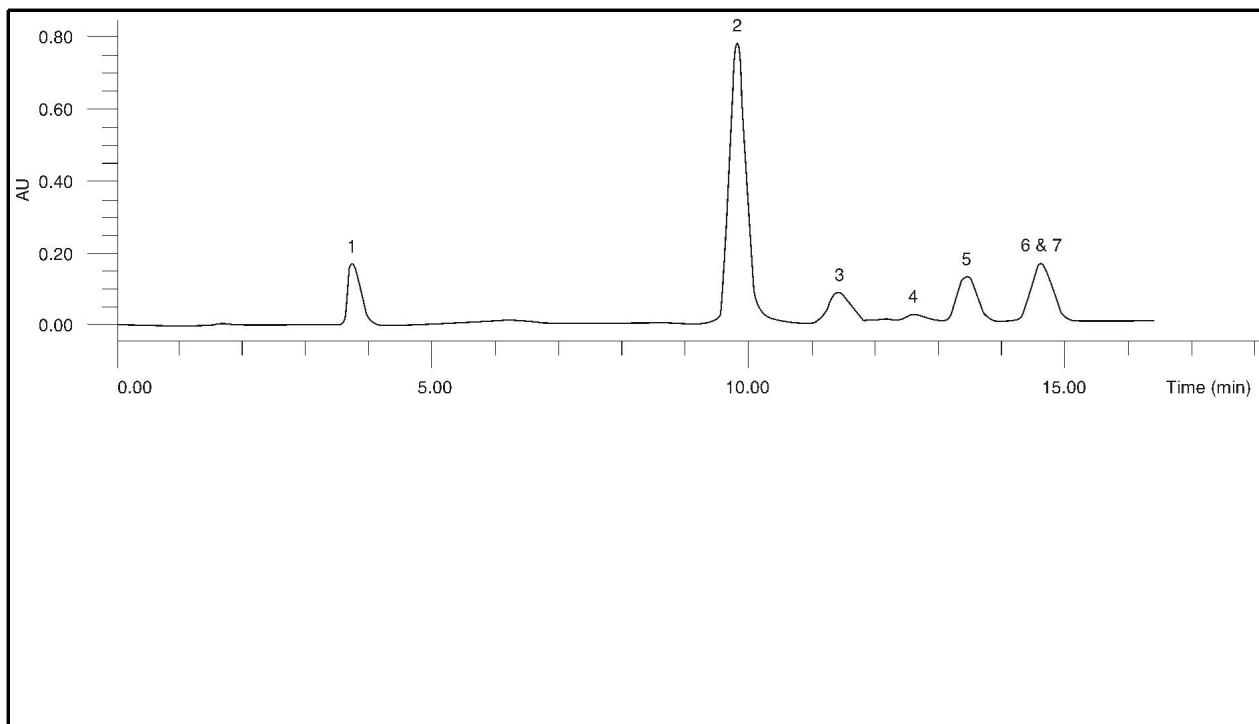

Column: Luna 5 μ m C18 / ODS Original

Dimensions: 250 x 4.6 mm ID

Order No:

Elution Type: Isocratic

Eluent A: Methanol

Pct. A: 80

Eluent B: Water

Pct. B: 20

Flow Rate: 1.5 mL/min

Temperature: 22 °C

Detection: UV-Vis Abs. - Variable Wave. @ 240 nm

SAMPLE ANALYTE

- 1 *Benzophenone3*
- 2 *DCHP*
- 3 *Parsol MCX*
- 4 *Homosalate 2*
- 5 *Octyl salicylate*
- 6 *Homosalate 1*
- 7 *Homosalate*

SPE:

For more information call your Phenomenex Representative at (310) 212-0555

Copyright 2000 Phenomenex Inc.

Ref. No.: 1239/ALG/100