

Phthalate Esters on Envirosep-PP (EPA Method 606)

Page 1 of 1

Column: EnviroSep-PP

Dimensions: 125 x 3.2 mm ID

Order No: 00E-3029-R0

Elution Type: Gradient

Eluent A: Water



Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	30	70
	2	4	30	70
	3	5	0	100

Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 254 nm (ambient)

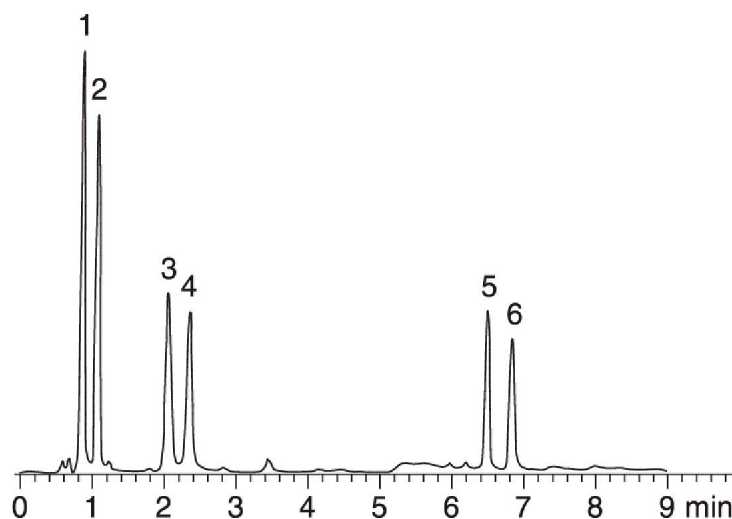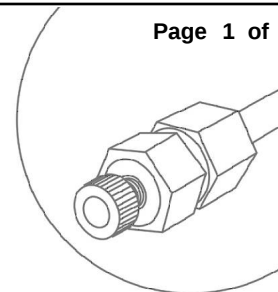

ANALYTES:

- 1 Dimethyl phthalate
- 2 Diethyl phthalate
- 3 Di-n-butyl phthalate
- 4 Butylbenzyl phthalate
- 5 Bis(2-ethylhexyl) phthalate
- 6 Di-n-octyl phthalate



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.

USA
tel.: (310) 212-0555
fax: (310) 328-7768
email: info@phenomenex.com

Puerto Rico
(800) 541-HPLC
(310) 328-7768
info@phenomenex.com

Canada
(800) 543-3681
(310) 328-7768
info@phenomenex.com

France
01 30 09 21 10
01 30 09 21 11
franceinfo@phenomenex.com

United Kingdom
01625-501367
01625-501796
ukinfo@phenomenex.com

Germany
06021-58830-0
06021-58830-11
anfrage@phenomenex.com

New Zealand
09-4780951
09-4780952
info@phenomenex.co.nz

Australia
02-9428-6444
02-9428-6445
info@phenomenex.com.au