



SEC of Proteins on Biosep-SEC-S3000

Column: BioSep-SEC-S 3000**Dimensions:** 600 x 7.8 mm ID**Order No:** 00K-2146-K0**Elution Type:** Isocratic**Eluent A:** 50 mM NaH₂PO₄ buffer, pH 6.8**Pct. A:** 100**Flow Rate:** 0.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs. - Variable Wave. @ 280 nm (ambient)

Page 1 of 1

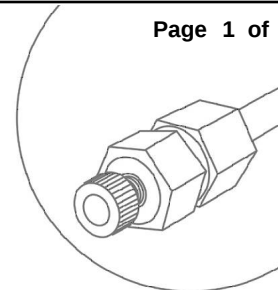

Products used in this Application:

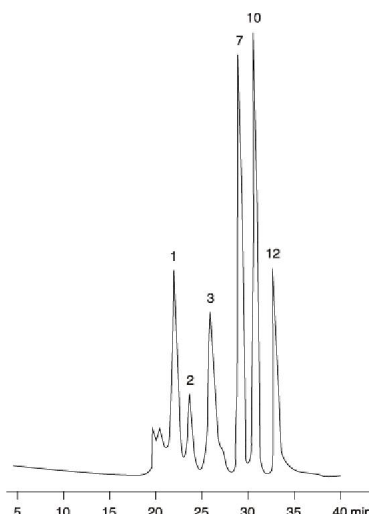Bio.Sep[™]

ANALYTES:

- 1 Thyroglobulin, 670K
- 2 IgA, 300K
- 3 IgG, 150K
- 4 LDH, 143K
- 5 Phosphorylase b, 97K
- 6 BSA, 68K
- 7 Ovalbumin, 44K
- 8 Horse radish peroxidase, 40K
- 9 Carbonic anhydrase, 30K
- 10 Trypsin inhibitor, 20.1K
- 11 Myoglobin, 17K
- 12 Lysozyme, 14.4K
- 13 Insulin, 5.7K
- 14 Cyanocobalamin, 1.35K

www.phenomenex.comPhenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.