

St. John's Wort on Synergi Max-RP

Page 1 of 1

Column: Synergi 4μ Max-RP 80Å

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4337-E0

Elution Type: Isocratic

Eluent A: 100mM Ammonium acetate

Pct. A: 20

Eluent B: Acetonitrile

Pct. B: 50


Eluent C: Methanol

Pct. C: 30

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs. - Variable Wave. @ 590 nm (ambient)

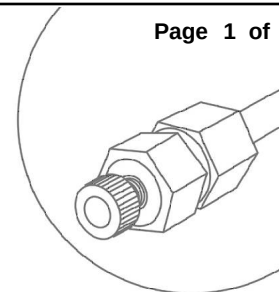

Products used in this Application:

synergi™

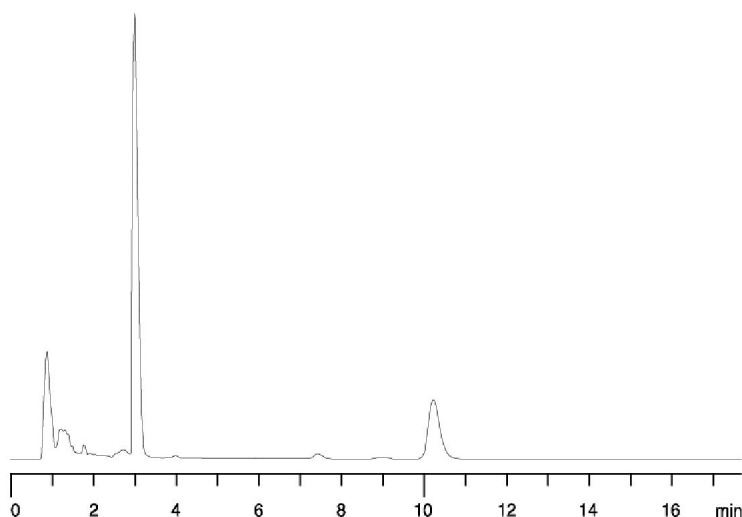

ANALYTES:

- 1 Pseudohypericin
- 2 Hypericin



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.