

Cocaine and Metabolites on Luna SCX

Column: Luna 5u SCX 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Isocratic

Eluent A: 50mM KH₂PO₄, pH 2.5

Pct. A: 65


Eluent B: Acetonitrile

Pct. B: 35

Flow Rate: 2 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

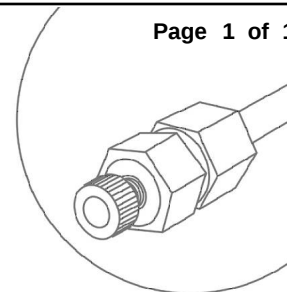

Page 1 of 1

Products used in this Application:

LUNA®

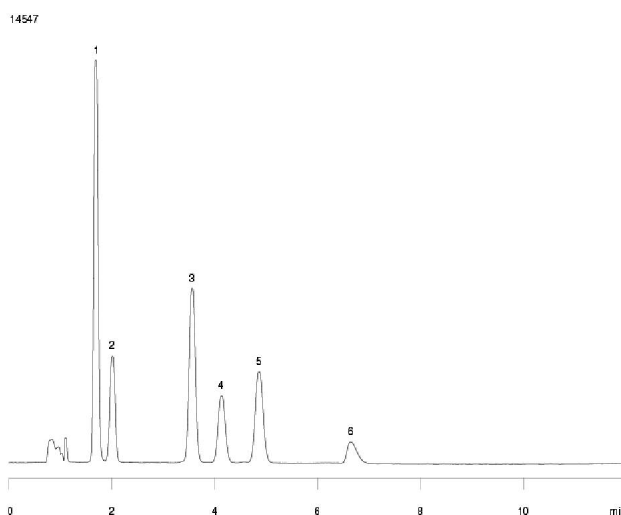

ANALYTES:

- 1 m-Hydroxybenzoyl Ecgonine
- 2 Norbenzoyl Ecgonine
- 3 Norcocaine
- 4 Benzoyl Ecgonine Isopropyl Ester
- 5 Cocaine
- 6 Ecgonine Methyl Ester

