

Oxycodone Derivatives on Luna SCX

Column: Luna 5u SCX 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Isocratic

Eluent A: 50mM KH₂PO₄, pH 2.5

Pct. A: 50

Eluent B: Acetonitrile


Pct. B: 50

Flow Rate: 2 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Page 1 of 1

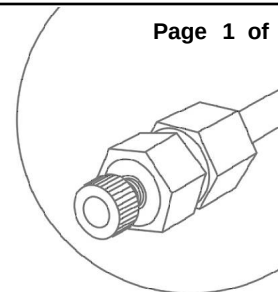

Products used in this Application:

LUNA®

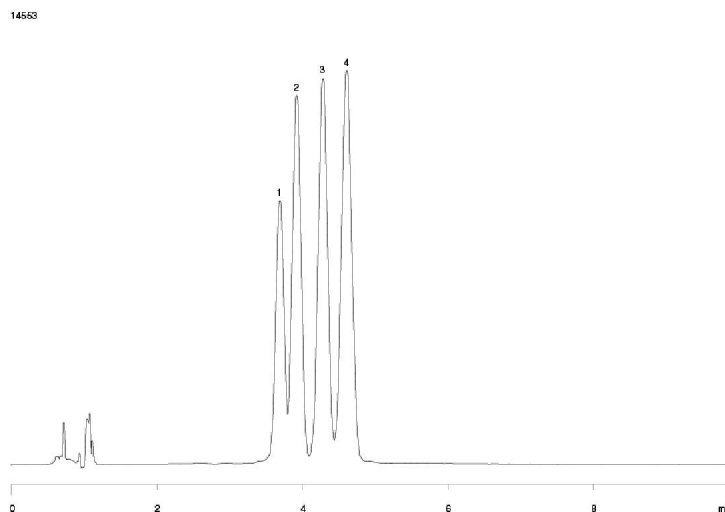

ANALYTES:

- 1 Noroxymorphone
- 2 Noroxycodone
- 3 Oxymorphone
- 4 Oxycodone

