

Children's Tylenol Cold on Luna SCX

Column: Luna 5u SCX 100A

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Isocratic

Eluent A: 50mM KH₂PO₄, pH 2.5

Pct. A: 35

Eluent B: Acetonitrile


Pct. B: 65

Flow Rate: 1.5 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs. - Variable Wave. @ 210 nm (ambient)

Page 1 of 1

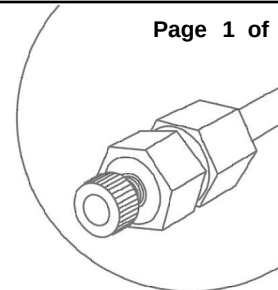

Products used in this Application:

LUNA®

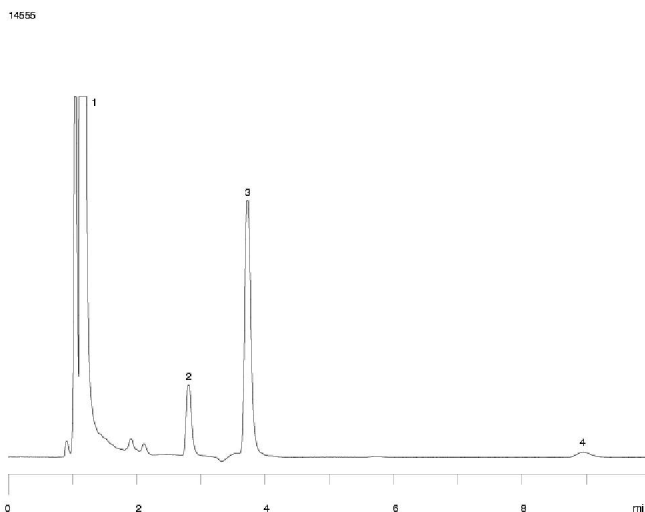

ANALYTES:

- 1 Acetaminophen
- 2 Pseudoephedrine
- 3 Dextromethorphan
- 4 Chlorpheniramine



www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country, contact Phenomenex USA, International Department by telephone, fax or e-mail: international@phenomenex.com.